

Cover Sheet for

**ENVIRONMENTAL CHEMISTRY METHOD**

***Pesticide Name:*** Glyphosate

***MRID #:*** 408816-01

***Matrix:*** Water

***Analysis:*** HPLC/FLD

This method is provided to you by the Environmental Protection Agency's (EPA) Environmental Chemistry Laboratory (ECL). This method *is not* an EPA method but one which was submitted to EPA by the pesticide manufacturer to support product registration. EPA recognizes that the methods may be of some utility to state, tribal, and local authorities, but makes no claim of validity by posting these methods. Although the Agency reviews *all* Environmental Chemistry Methods submitted in support of pesticide registration, the ECL evaluates only about 30% of the currently available methods. Most methods perform satisfactorily but some, particularly the older methods, have deficiencies. Moreover, the print quality of the methods varies considerably because the methods originate from different sources. Therefore, the methods offered represent the best available copies.

If you have difficulties in downloading the method, or further questions concerning the methods, you may contact Elizabeth Flynt at 228-688-2410 or via e-mail at [flynt.elizabeth@epa.gov](mailto:flynt.elizabeth@epa.gov).



**Analytical Residue Method for N-Phosphonomethyl Glycine and  
Aminomethylphosphonic Acid in Environmental Water**

Method No. 86-63-1

**Scope**

The procedures given can determine trace quantities of N-(phosphonomethyl) glycine (glyphosate) and aminomethylphosphonic acid (AMPA) in parts per billion range in environmental water by HPLC.

**Summary**

Environmental water can be analyzed for glyphosate and AMPA by concentration and injection into a high pressure liquid chromatograph equipped with an o-phthalaldehyde (OPA) post-column reactor (PCR) and a fluorescence detector.

**Sensitivity**

The sensitivity of the method is 0.5 ppb after concentration of 250 mL of environmental water for both glyphosate and AMPA. If the environmental water is injected straight into the HPLC-PCR system, a concentration of 25 ppb can be determined directly by external standard calibration.

**Apparatus and Equipment**

High vacuum pump, Sargent-Weich Model 1400B or equivalent

Rotary Evaporator

Cold finger condenser (45 cm long) filled with dry ice

3 mL disposable syringes with Luer-Lok® tips (Becton-Dickinson)

Arcodisc® disposable filter assembly, 0.45 µg pore size

6 inch cotton tipped applicators

500 mL round bottom flask

Vials for autosampler

## Reagents

200

## A. Chemicals and Materials

Deionized water

Hydrochloric acid, reagent grade

Potassium dihydrogen phosphate (HPLC grade)

Methanol (HPLC grade)

Phosphoric acid, concentrated (HPLC grade)

Disodium ethylenediamine tetraacetate dihydrate, certified A.C.S.

## B. Fortification Solutions

## 1. Glyphosate and AMPA

Weigh and dissolve 0.1000 g of N-(phosphonomethyl) glycine (Glyphosate) and 0.1000 g of aminomethylphosphonic acid (AMPA) in 1000 mL of deionized water. This concentrate contains 100 micrograms of glyphosate and AMPA per milliliter.

Standard solutions for fortifications are made by diluting the standard concentrate with deionized water as follows:

## From 100 Microgram Concentrate

Milliliters Concentrate	Standard Dilution	Concentration Micrograms/mL
50.0	100.0	50.0
30.0	100.0	30.0
20.0	100.0	20.0
10.0	100.0	10.0
5.00	100.0	5.0
3.00	100.0	3.0
2.00	100.0	2.0
1.00	100.0	1.0
0.50	100.0	.50
0.25	100.0	.25

## 2. Source of Standard Material Used in Fortifications and HPLC Standard Solutions.

Glyphosate reference N-(phosphonomethyl) glycine > 99

AMPA reference (aminomethylphosphonic acid), 99%

### C. HPLC Standards

HPLC standard solutions are prepared in a similar manner as above except that 0.001 M disodium EDTA solution is used in place of deionized water to make dilutions. The 100 and 10.0 microgram per milliliter stock solutions are made as for the fortification solutions with 0.001 M disodium EDTA solution. The solutions for HPLC standards are made by diluting the 10.0 microgram per milliliter stock with 0.001 M disodium EDTA solution as follows:

#### From 10.0 Microgram Dilution

Milliliters Concentrate	Standard Dilution	Concentration Micrograms/mL
10.0	100.0	1.00
5.00	100.0	0.50
1.00	100.0	0.10
0.50	100.0	.050
0.25	100.0	.025

### Procedure

#### Sample Preparation

Thaw frozen water sample, shake thoroughly, and transfer approximately half of the 250 mL sample (filter through glass wool or filter paper if there is a lot of suspended particles in the sample) into a 500 mL round bottom flask. For recovery samples, analyte fortifications are made at this stage. Add 5 mL of concentrated hydrochloric acid to the sample in the flask and 5 mL to the sample remaining in the bottle. Concentrate the sample on a rotary-film evaporator by slowly increasing the temperature of the water bath from 20° C to 60° C. Before the sample is evaporated to dryness, add the remainder of the sample and rinse the bottle twice with approximately 5 mL of deionized water, adding the rinses to the flask also. Concentrate the sample to dryness and remove the final traces of moisture with a stream of dry nitrogen if necessary. Combine the residue with 2.9 mL of HPLC buffer 0.005 M  $KH_2PO_4$  in 4% methanol/deionized water adjusted to pH 2.1 with concentrated phosphoric acid) using a long cotton swab

to remove any solid residues from the bottom and a pipette to rinse the sides of the flask. Leaving the swab in the flask, add 0.1 mL of 0.03 M disodium EDTA solution and mix thoroughly by rinsing the sides. Transfer the solution into a disposable syringe and filter through an attached 0.45  $\mu\text{m}$  size pore membrane filter. The sample is now ready for quantification of glyphosate and AMPA by using the HPLC - OPA post-column reactor system.

Samples that need to be diluted to remain within the standard concentration range must be diluted with a 0.001 M EDTA solution.

**Note:** Standards used for quantification of concentrated samples containing EDTA must also contain EDTA. If environmental water is injected straight into HPLC-PCR system, filter through a 0.45  $\mu\text{m}$  filter and compare to standards dissolved in deionized water.

### HPLC OPA Post-Column Reactor System

Glyphosate and aminomethylphosphonic acid (AMPA) may be separated and detected utilizing a high pressure liquid chromatograph and a post column reaction specific for primary amines. Glyphosate is oxidized with calcium hypochlorite and the product (glycine) and AMPA are coupled with O-phthalaldehyde in the presence of mercaptoethanol (OPA-MERC) to give fluorophors detected by a fluorometer with excitation at 340 nm and emission measured at 455 nm. The components needed for the construction of this detection system are outlined below. A general schematic is presented in Figure 1 and several general comments on the assembly are included.

Comments on the assembly are included. The system is designed for the detection of glyphosate and AMPA. The reaction conditions are optimized for the detection of these compounds. The system is designed for the detection of glyphosate and AMPA. The reaction conditions are optimized for the detection of these compounds. The system is designed for the detection of glyphosate and AMPA. The reaction conditions are optimized for the detection of these compounds.

**Equipment and Supplies**

- Solvent Buffer Reservoir
- HPLC pumps (two required)
- Temperature controller for HPLC column and oxidation reaction coil (Kratos URA 200 or equivalent)
- $Ca(OCl)_2$  solution reservoir
- OPA solution reservoir
- Post Column Derivatization System (Kratos URS 051 dual pump Reagent Delivery Module or equivalent)
- 1 mL reaction coil (2)
- stainless steel tubing  $\left\{ \begin{array}{l} 1/16 \text{ inch} \times \text{O.D.} 0.020 \text{ inch LD.} \\ 1/16 \text{ inch} \times \text{O.D.} 0.010 \text{ inch LD.} \end{array} \right.$
- Fluorescence Spectrometer (Perkin-Elmer LS-4 or equivalent)
- Solvent filter (Millipore GSWP 04700 or equivalent)
- Electronic integrator or computer

**Reagents**

- Deionized water
- Potassium dihydrogen phosphate (HPLC grade)
- Methanol (HPLC grade)
- Phosphoric acid, concentrated (HPLC grade)
- Fluoraldehyde<sup>®</sup>, o-phthalaldehyde reagent solution (Pierce Chemical Company)
- Disodium ethylenediamine tetraacetate dihydrate, certified A.C.S. grade
- Alternate Solution for Fluoraldehyde<sup>®</sup>
  1. Boric acid
  2. 2-Mercaptoethanol
  3. Potassium hydroxide, 45% solution
  4. Brig<sup>®</sup>35, 30% solution
  5. Deionized water, with resistance of 18 megohm/cm
  6. Fluoropa o-phtha.a.dehyde, Dionex) Pierce Chemical Company
  7. Methanol

Calcium hypochlorite, certified, 70.89% available chlorine,

Fisher Scientific Company

Sodium Chloride, analytical reagent

Sodium hydroxide, certified A.C.S. reagent

### Solution Preparation

#### 1. HPLC Buffer

Prepare 0.005 M potassium dihydrogen phosphate ( $KH_2PO_4$ ) by dissolving 2.72 g in four liters of 4% methanol/deionized water. This solution is adjusted to pH 2.1 with concentrated phosphoric acid. Normal HPLC degassing procedures are followed, and the solution is filtered through a 0.22  $\mu$ m filter.

#### 2. Oxidative Solution

In a 500 mL volumetric flask, dissolve 0.5 g  $Ca(OCl)_2$  in 500 mL of deionized water using a magnetic stirrer at high speed for 45 minutes. In a 1 liter volumetric flask, dissolve 1.36 g  $KH_2PO_4$ , 11.6 g  $NaCl$ , and 0.4 g  $NaOH$  in 500 mL solution) and dilute to 1 L with deionized water, mixing well. Filter the solution through a 0.22  $\mu$ m filter.

#### 3. EDTA Solutions

##### a. 0.03 M EDTA Solution

Dissolve 11.17 g of disodium ethylenediamine tetraacetate dihydrate in 1 L deionized water using a magnetic stirrer and filter through a 0.22  $\mu$ m filter.

##### b. 0.001 M EDTA Solution

Dissolve 0.37 g of disodium ethylenediamine tetraacetate dihydrate in 1 L deionized water using a magnetic stirrer and filter through a 0.22  $\mu$ m filter.

#### 4. Alternate OPA Solution

a. Dissolve 25 g boric acid in 950 mL deionized water, using a magnetic stirrer.

b. While monitoring the pH with a pH meter, titrate with the potassium hydroxide solution (approximately 30 mL will be required) to a final pH of 10.40  $\pm$  0.2.

c. Filter through a 0.22  $\mu$ m filter.

d. Add 3 mL 30% Brig 35 solution.

e. Add 2.0 mL 2-mercaptoethanol.



- f. Dissolve 800 mg Fluoropa in 10 mL methanol at room temperature using gentle swirling.
- g. Add the methanol-Fluoropa solution.

#### Storage Stability

Fluoropa and 2-mercaptoethanol are subject to atmospheric oxidation and these oxidation products can contribute to increased background fluorescence. Thus, unless the reagent solution is protected from atmospheric oxygen, it should be prepared fresh daily. The solution can be stored in closed glass bottles under atmospheric conditions at 4° C for up to two weeks without appreciable increases in background fluorescence or it can be stored under nitrogen for indefinite periods.

*Note: The Fluoraldehyde® (Pierce OPA reagent solution) is a specially formulated OPA solution which has outstanding shelf life with no increase in background fluorescence with time.*

#### HPLC Conditions

Pre-Column: RP-18 Spheri-10, 3.6 cm x 4.6 mm I.D., guard column, Brownlee Labs Inc., Santa Clara, California.

Column: Aminex A-9, 30 cm x 4.6 mm I.D., Bio-Rad Laboratories, Richmond, California.

Column Temperature: 50° C

Buffer Flow: 0.5 mL/min

Oxidation solution flow rate: 0.5 mL/min

OPA solution flow rate: 0.5 mL/min

Pressure: approximately 1500 psi Buffer

Fluorometer settings: { Excitation 340 nm  
Emission 455 nm

Injection volume: 200 µL

See Figure 1 for an outline of the system.

### General Comments

The system is designed to permit the continuous analysis of samples. A pre-column or guard column is installed to protect the analytical column. This pre-column can be changed periodically or when necessary.

The analytical column should be packed using the HPLC buffer as eluant. The column end fittings should accommodate 1/16 x 0.001 inch I.D. tubing and must be capable of withstanding a backpressure of 4000 psi.

When starting, turn all pumps and the detector on for thirty minutes prior to use. If an air bubble becomes trapped in the detector cell, it can be removed by disconnecting the waste line from the detector and alternately drawing and forcing liquid through the cell with a syringe containing water or methanol, until the bubble is removed. When shutting the system down, it is advisable to flush the system (pumps and fluorometer) with deionized water.

The relative peak response of glyphosate to AMPA can be adjusted by optimizing the flow rates of hypochlorite and OPA. Increasing the hypochlorite flow rate will decrease the AMPA response relative to the glyphosate response.

A buffer of pH 2.1 was used for analyzing glyphosate and aminomethylphosphonic acid in water samples. Sometimes a close peak interfering with the integration can be further separated by changing the pH of the buffer. Should it become necessary to change the pH of the solvent buffer to effect a better separation, the hypochlorite and OPA flows should also be adjusted to optimize the greatest response for glyphosate and AMPA. The pH and the amount of excess hypochlorite and OPA solutions present do effect the reactions that produce the fluorogenic responses.

The  $Ca(OCl)_2$  reaction coil can be run at room temperature but temperature should be constant throughout the run. Variation of temperature by several degrees can effect the response of the glyphosate peak. To prevent problems from temperature variation, the  $Ca(OCl)_2$  reaction coil should be placed in a temperature controlled module such as the Kratos URA 100 or 200 reaction chamber/column temperature control module. A temperature of 40-C gives the optimum response for both glyphosate and AMPA.

Using the HPLC conditions given above, the analysis time between injections is approximately 30 minutes, depending on the quality of the column.

**Quantitation and Quality Control**

Sample quantitation is based on the relative peak height or peak area of the sample to standard peak heights or areas across the range of expected sample concentration. Laboratory fortification of field samples serves to estimate the recovery of the analytical method. A minimum of two laboratory fortifications are run with each sample set. Upon completion of all samples the recovery data from each set is compiled and averaged to generate a recovery estimate for the method over the entire protocol.

**Method Validation**

Appropriate method validation in conjunction with quality control is necessary to insure the integrity of the data generated.

1943

Presented to the Board of Directors

1943

The following is a summary of the financial results of the company for the year ended December 31, 1943. The company has achieved a net income of \$100,000, which is a significant increase over the previous year. This is due to the successful completion of several major projects and the implementation of cost-saving measures. The company's assets have increased by \$200,000, and its liabilities have decreased by \$50,000. The overall financial position of the company is strong and stable.

Respectfully,  
[Signature]

The Board of Directors has reviewed the financial statements and has approved the same. The company's performance has been excellent, and the Board is confident in the company's future prospects. The Board has also approved the dividend payment of \$50,000 to the shareholders.

1943

MSL-8332

205

EXAMPLE CHROMATOGRAMS

SECRET

001

SECRET


SECRET

MSL 8332


Representative Chromatograms of Flowing Stream  
HPLC Analysis - April 11, 1988

1, MI


211

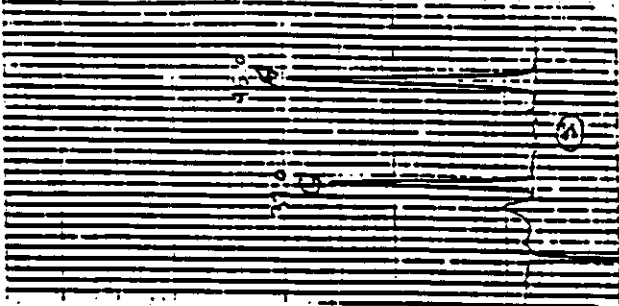


Glyphosate standards  
0.50  $\mu$ g each

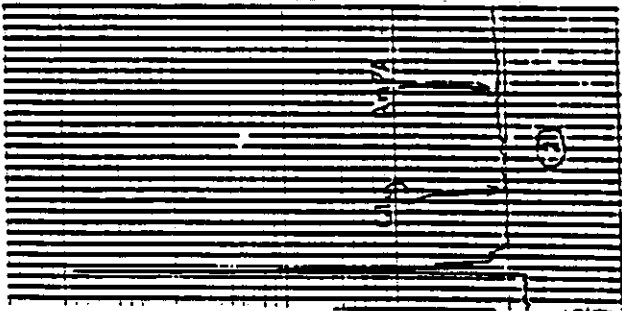

Control Stream Water, -1 DAT  
<0.001 ppm Glyphosate, <0.001 ppm AMPA

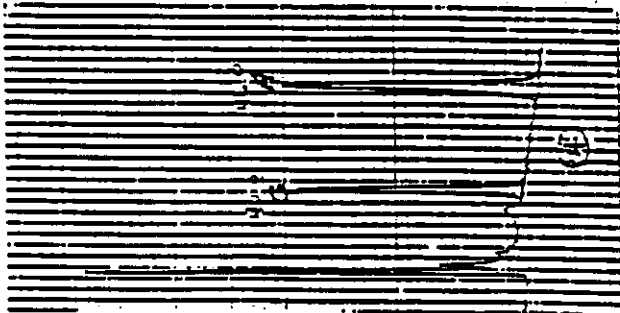

Control Stream Water, -1 DAT  
Fortified @ 0.010 ppm each  
Recoveries: 107.71% Glyphosate  
101.66% AMPA

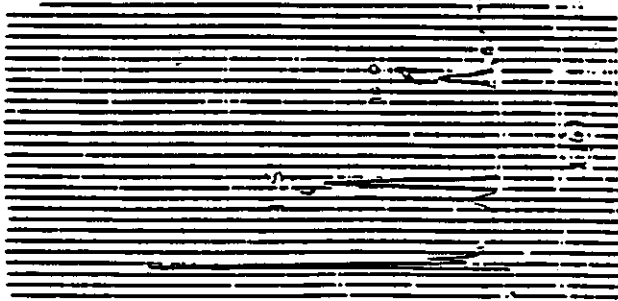


Treated Stream Water, 3 DAT  
0.008 ppm Glyphosate, 0.002 ppm AMPA

0.001


MSL 8332

Representative Chromatograms of Pond Sediment from Clarence, MO  
HPLC Analysis - December 13, 1987


213

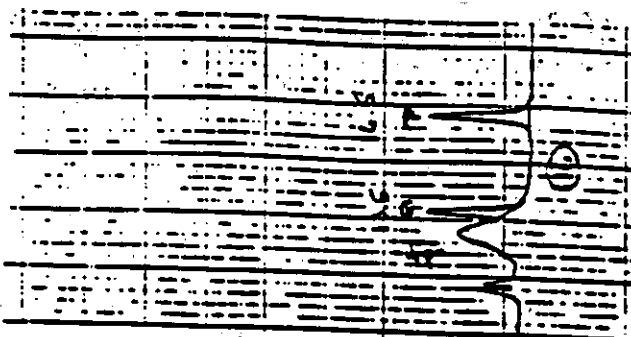


Glyphosate and AMPA Standards  
0.5 µg/mL each

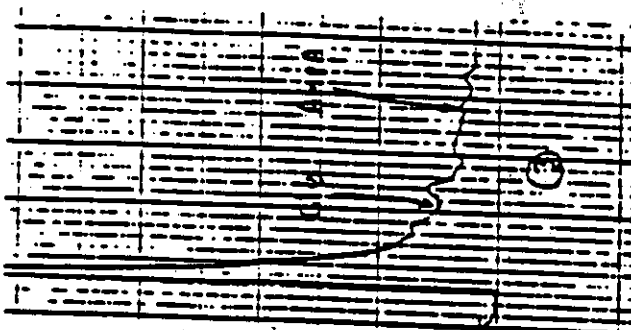

Control Sediment, Missouri, 7 DAT  
<0.05 ppm glyphosate, <0.05 ppm AMPA

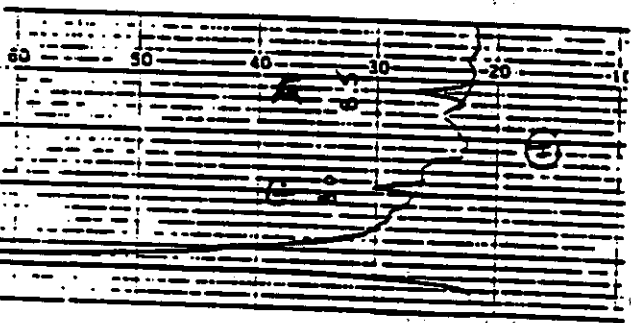

Control Sediment, Missouri, 7 DAT  
Fortified @ 0.05 ppm each  
recoveries: 37% Glyphosate, 80% AMPA

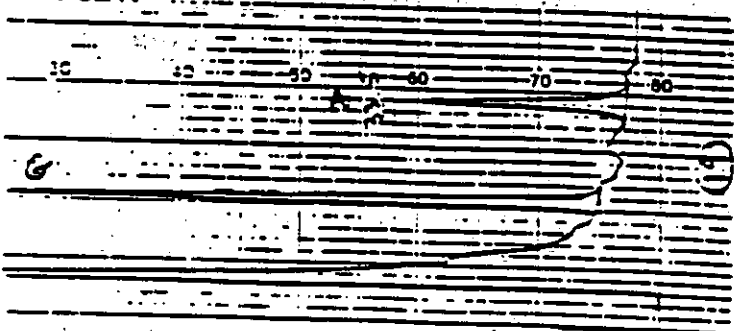

Treated Sediment, Missouri, 7 DAT  
0.79 ppm glyphosate, 0.22 ppm AMPA

MSL 8332  
HPLC Analysis - December 13, 1987  
Peak from previous collection




MSL #8332

Representative Chromatograms of Pond Sediment from Corvallis, OR  
HPLC Analysis - May 9, 1988


213




Glyphosate and AMPA Standards  
0.50  $\mu\text{g/mL}$  each

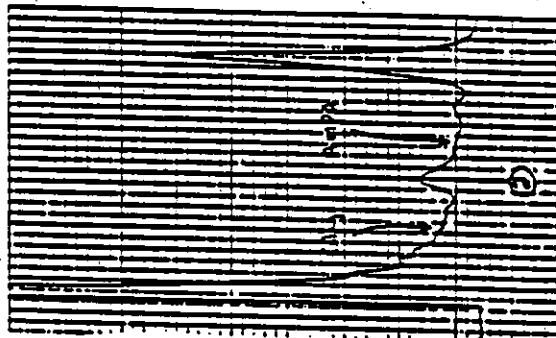

Control Pond Sediment, -9 DAT  
<0.05 ppm Glyphosate, <0.05 ppm AMPA

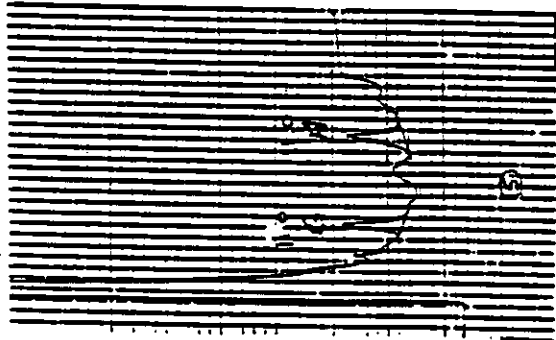

Control Pond Sediment, -9 DAT  
Fortified @ 0.10 ppm each  
Recoveries: 72.33% Glyphosate  
70.67% AMPA

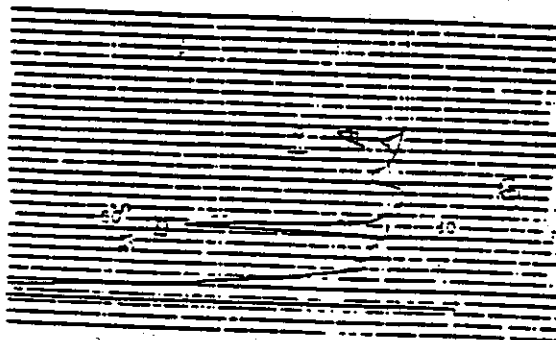

Treated Pond Sediment, 14 DAT  
0.38 ppm Glyphosate, 0.09 ppm AMPA


MSL 8332

Representative Chromatograms of Flowing Stream Sediment from Cuckbert, GA  
HPLC Analysis - May 5, 1988


214




Glyphosate and AMPA Standards  
0.25  $\mu\text{g}/\text{mL}$  each

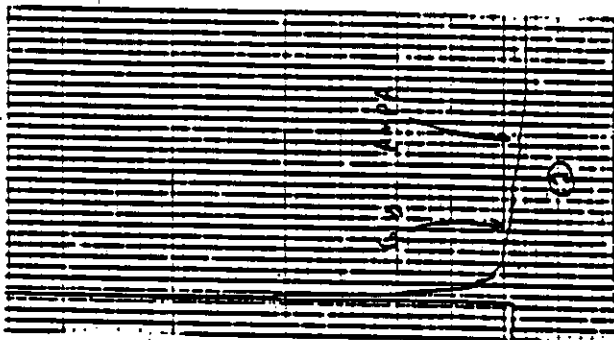

Control Stream Sediment, -1 DAT  
<0.05 ppm Glyphosate, <0.05 ppm AMPA

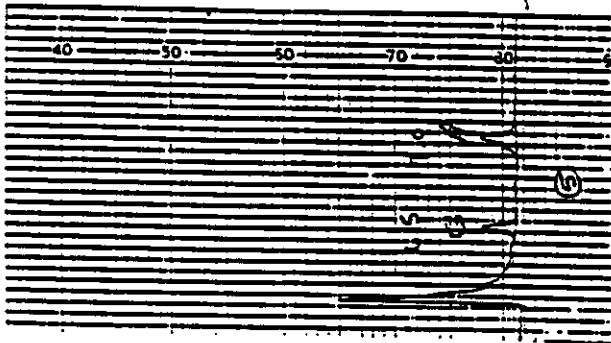

Control Stream Sediment, -1 DAT  
Fortified @ 0.05 ppm each  
Recoveries: 95.08% Glyphosate  
88.38% AMPA

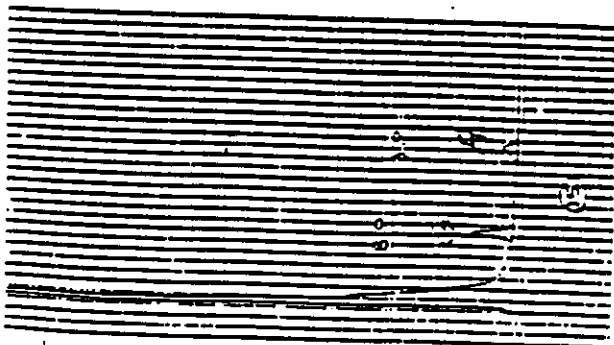

Treated Stream Sediment, 58 DAT  
0.08 ppm Glyphosate, <0.05 ppm AMPA