

Cover Sheet for

**ENVIRONMENTAL CHEMISTRY METHOD**

***Pesticide Name:*** Milbemectin

***MRID # :*** 454147-04

***Matrix:*** Soil

***Analysis:*** HPLC/UV

This method is provided to you by the Environmental Protection Agency's (EPA) Environmental Chemistry Laboratory (ECL). This method is *not* an EPA method but one which was submitted to EPA by the pesticide manufacturer to support product registration. EPA recognizes that the methods may be of some utility to state, tribal, and local authorities, but makes no claim of validity by posting these methods. Although the Agency reviews *all* Environmental Chemistry Methods submitted in support of pesticide registration, the ECL evaluates only about 30% of the currently available methods. Most methods perform satisfactorily but some, particularly the older methods, have deficiencies. Moreover, the print quality of the methods varies considerably because the methods originate from different sources. Therefore, the methods offered represent the best available copies.

If you have difficulties in downloading the method, or further questions concerning the methods, you may contact Elizabeth Flynt at 228-688-2410 or via e-mail at [flynt.elizabeth@epa.gov](mailto:flynt.elizabeth@epa.gov).

VALIDATION OF ANALYTICAL METHODS FOR THE DETERMINATION  
OF MILBEMECTIN AND ITS POTENTIAL DEGRADATION PRODUCTS  
IN FIVE SOIL TYPES

WILDLIFE INTERNATIONAL LTD. PROJECT NUMBER: 420C-106

AUTHORS:

Timothy Z. Kendall, M.S.  
Willard B. Nixon, Ph.D.

SUBMITTED TO:

Sankyo Co., Ltd  
Agrochemicals Division  
No. 7-12, Ginza 2-Chome  
Chuo-ku, Tokyo 104, Japan

*Wildlife International, Ltd.*

8598 Commerce Drive  
Easton, MD 21601  
(410) 822-8600

Total Pages. 65

WILDLIFE INTERNATIONAL LTD

PROJECT NO.: 420C-106

- 2 -

## TABLE OF CONTENTS

TITLE/COVER .....	1
TABLE OF CONTENTS .....	2
1.0 INTRODUCTION .....	6
1.1 Scope.....	6
1.2 Principle .....	6
2.0 APPARATUS.....	8
3.0 <u>REAGENTS</u> .....	8
4.0 REFERENCE STANDARDS IDENTIFICATION.....	9
5.0 ANALYTICAL STANDARD SOLUTION PREPARATION .....	10
5.1 Stock Solution Preparation - M.A <sub>3</sub> and M.A <sub>4</sub> .....	10
5.2 Stock Solution Preparation - 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	10
5.3 Stock Solution Preparation - 27-hydroxy - M.A <sub>3</sub> and 27-hydroxy -M.A <sub>4</sub> .....	10
5.4 Stock Solution Preparation - 27-keto - M.A <sub>3</sub> and 27-keto -M.A <sub>4</sub> .....	10
5.5 Combined Stock Solution - M.A <sub>3</sub> , M.A <sub>4</sub> , 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	11
5.6 Calibration Standard Solutions of Anhydro- M.A <sub>3</sub> and Anhydro-M A <sub>4</sub> .....	11
5.7 Calibration Standard Solutions of Anhydro-27-Hydroxy-M.A <sub>3</sub> and Anhydro-27-Hydroxy-M.A <sub>4</sub> .....	11
5.8 Calibration Standard Solutions of 27-Keto-M.A <sub>3</sub> and 27-Keto-M.A <sub>4</sub> .....	12
6.0 ANALYTICAL METHOD .....	12
6.1 Sample Preparation .....	12
6.2 Laboratory Fortification.....	12
6.3 Sample Extraction .....	13
6.4 Solid-Phase Extraction.....	14
6.5 Separation of Milbemectin and its Potential Degradation Products.....	14
6.6 Derivatization .....	15
6.7 Separation of M.A <sub>3</sub> /M.A <sub>4</sub> from 8,9Z-M.A <sub>3</sub> /8,9Z-M.A <sub>4</sub> .....	15
6.8 Fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy-M.A <sub>3</sub> /27-hydroxy-M.A <sub>4</sub> and 27-keto-M.A <sub>3</sub> /27-keto-M.A <sub>4</sub> .....	16
6.9 HPLC Determination.....	17
7.0 CHROMATOGRAPHIC CONDITIONS .....	17
7.1 HPLC Conditions for the Separation of Anhydro-M A <sub>3</sub> and Anhydro-M.A <sub>4</sub> Derivatives.....	17
7.2 HPLC Conditions for Fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> and 8,9Z-M.A <sub>3</sub> /8,9Z-M.A <sub>4</sub> .....	17
7.3 Isocratic HPLC Operational Parameters For The Analysis of Anhydro-M.A <sub>3</sub> and Anhydro-M A <sub>4</sub> .....	18
7.4 Normal Phase HPLC conditions for the Analysis of 27 Keto-M.A <sub>3</sub> and 27 Keto-M.A <sub>4</sub> .....	18

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 3 -

8.0	CALIBRATION .....	19
9.0	CALCULATIONS .....	19

## TABLES

Table 1.	Typical Gradient HPLC Conditions for the Analysis of M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy-M.A <sub>3</sub> /27-Hydroxy-M.A <sub>4</sub> .....	21
Table 2.	Typical Normal Phase HPLC Conditions for the Fractionation of 8,9Z-M.A <sub>3</sub> , 8-9Z-M.A <sub>4</sub> , M.A <sub>3</sub> and M.A <sub>4</sub> .....	22
Table 3.	Typical Normal Phase HPLC Conditions for the Analysis of 27-Keto M.A <sub>3</sub> , and 27-Keto M.A <sub>4</sub> .....	23
Table 4.	Method Validation Recoveries of Fractionated M.A <sub>3</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	24
Table 5.	Method Validation Recoveries of Fractionated M.A <sub>4</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	25
Table 6.	Method Validation Recoveries of Fractionated 8,9Z-M.A <sub>3</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	26
Table 7.	Method Validation Recoveries of Fractionated 8,9Z-M.A <sub>4</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 8,9Z-M.A <sub>3</sub> and 8,9Z-M.A <sub>4</sub> .....	27
Table 8.	Method Validation Recoveries of M.A <sub>3</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy M.A <sub>3</sub> , 27-Hydroxy M.A <sub>4</sub> , 27-Keto M.A <sub>3</sub> , 27-Keto M.A <sub>4</sub> .....	28
Table 9.	Method Validation Recoveries of M.A <sub>4</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy M.A <sub>3</sub> , 27-Hydroxy M.A <sub>4</sub> , 27-Keto M.A <sub>3</sub> , 27-Keto M.A <sub>4</sub> .....	29
Table 10.	Method Validation Recoveries of 27-Hydroxy M.A <sub>3</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy M.A <sub>3</sub> , 27-Hydroxy M.A <sub>4</sub> , 27-Keto M.A <sub>3</sub> , 27-Keto M.A <sub>4</sub> .....	30
Table 11.	Method Validation Recoveries of 27-Hydroxy M.A <sub>4</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy M.A <sub>3</sub> , 27-Hydroxy M.A <sub>4</sub> , 27-Keto M.A <sub>3</sub> , 27-Keto M.A <sub>4</sub> .....	31
Table 12.	Method Validation Recoveries of 27-Keto M.A <sub>3</sub> from New York Sandy Loam Fortified with M.A <sub>3</sub> , M.A <sub>4</sub> , 27-Hydroxy M.A <sub>3</sub> , 27-Hydroxy M.A <sub>4</sub> , 27-Keto M.A <sub>3</sub> , 27-Keto M.A <sub>4</sub> .....	32

Table 13. Method Validation Recoveries of 27-Keto M.A. <sub>4</sub> from New York Sandy Loam Fortified with M.A. <sub>3</sub> , M.A. <sub>4</sub> , 27-Hydroxy M.A. <sub>3</sub> , 27-Hydroxy M.A. <sub>4</sub> , 27-Keto M.A. <sub>3</sub> , 27-Keto M.A. <sub>4</sub> .....	33
---	----

## FIGURES

Figure 1. Analytical method flow chart for the analysis of Milbemectin and its metabolites in soil.....	34
Figure 2. Analytical method flow chart for the fractionation of the parent and photoproducts of milbemectin.....	36
Figure 3. Flow chart for Processing of Milbemectin and Potential Metabolites in Soil.....	37
Figure 4. A representative calibration curve for M.A. <sub>3</sub> .....	38
Figure 5. A representative calibration curve for M.A. <sub>4</sub> .....	39
Figure 6. A representative calibration curve for 27-Hydroxy M.A. <sub>3</sub> .....	40
Figure 7. A representative calibration curve for 27-Hydroxy M.A. <sub>4</sub> .....	41
Figure 8. A representative calibration curve for 27-Keto M.A. <sub>3</sub> .....	42
Figure 9. A representative calibration curve for 27-Keto M.A. <sub>4</sub> .....	43
Figure 10. A representative chromatogram of a 5.00 µg a.i./L 27-Keto M.A. <sub>3</sub> and 27-Keto M.A. <sub>4</sub> combined standard.....	44
Figure 11. A representative chromatogram of a 50.0 µg a.i./L 27-Keto M.A. <sub>3</sub> and 27-Keto M.A. <sub>4</sub> combined standard.....	45
Figure 12. A representative chromatogram of a 0.50 µg a.i./L M.A. <sub>3</sub> , 27-Hydroxy M.A. <sub>3</sub> , M.A. <sub>4</sub> and 27-Hydroxy M.A. <sub>4</sub> combined standard.....	46
Figure 13. A representative chromatogram of a 5.00 µg a.i./L M.A. <sub>3</sub> , 27-Hydroxy M.A. <sub>3</sub> , M.A. <sub>4</sub> and 27-Hydroxy M.A. <sub>4</sub> combined standard.....	47
Figure 14. A representative normal phase chromatogram of the fractionation of 8,9Z-M.A. <sub>3</sub> and 8,9Z-M.A. <sub>4</sub> from M.A. <sub>3</sub> and M.A. <sub>4</sub> .....	48
Figure 15. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> analysis of the reagent blank from the fractionation of 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> from New York sandy loam.....	49
Figure 16. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> analysis of the matrix blank from the fractionation of 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> from New York sandy loam.....	50
Figure 17. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> fractionation analysis of New York sandy loam fortified at 2.00 µg a.i./L with 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> .....	51
Figure 18. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> fractionation analysis of New York sandy loam fortified at 10.0 µg a.i./L with 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> .....	52
Figure 19. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> fractionation analysis of New York sandy loam fortified at 50.0 µg a.i./L with 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> .....	53
Figure 20. A representative chromatogram of the M.A. <sub>3</sub> /M.A. <sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with 8,9Z-M.A. <sub>3</sub> /8,9Z-M.A. <sub>4</sub> and M.A. <sub>3</sub> /M.A. <sub>4</sub> .....	54

Figure 21.	A representative chromatogram of the M.A <sub>3</sub> /M.A <sub>4</sub> and 27-Hydroxy M.A <sub>3</sub> /27-Hydroxy M.A <sub>4</sub> analysis of the reagent blank from the fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> from New York sandy loam .....	55
Figure 22	A representative chromatogram of the M.A <sub>3</sub> /M.A <sub>4</sub> and 27-Hydroxy M.A <sub>3</sub> /27-Hydroxy M.A <sub>4</sub> analysis of the matrix blank from the fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> from New York sandy loam .....	56
Figure 23	A representative chromatogram of the M.A <sub>3</sub> /M.A <sub>4</sub> and 27-Hydroxy M.A <sub>3</sub> /27-Hydroxy M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 2.00 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> from New York sandy loam .....	57
Figure 24	A representative chromatogram of the M.A <sub>3</sub> /M.A <sub>4</sub> and 27-Hydroxy M.A <sub>3</sub> /27-Hydroxy M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 10.0 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> .....	58
Figure 25.	A representative chromatogram of the M.A <sub>3</sub> /M.A <sub>4</sub> and 27-Hydroxy M.A <sub>3</sub> /27-Hydroxy M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> .....	59
Figure 26.	A representative chromatogram of the 27-keto M.A <sub>3</sub> / 27-keto M.A <sub>4</sub> analysis of the reagent blank from the fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> from New York sandy loam .....	60
Figure 27.	A representative chromatogram of the 27-keto M.A <sub>3</sub> / 27-keto M.A <sub>4</sub> analysis of the reagent blank from the fractionation of M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> from New York sandy loam .....	61
Figure 28.	A representative chromatogram of the 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 5.00 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> .....	62
Figure 29	A representative chromatogram of the 27-keto M.A <sub>3</sub> / 27-keto M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 10.0 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> .....	63
Figure 30.	A representative chromatogram of the 27-keto M.A <sub>3</sub> / 27-keto M.A <sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with M.A <sub>3</sub> /M.A <sub>4</sub> , 27-hydroxy M.A <sub>3</sub> /27-hydroxy M.A <sub>4</sub> and 27-keto M.A <sub>3</sub> /27-keto M.A <sub>4</sub> .....	64

**APPENDIX**

Appendix I	Soil Characterization for Milbemectin.....	65
------------	--	----

WILDLIFE INTERNATIONAL LTD

PROJECT NO.: 420C-106

- 6 -

## 1.0 INTRODUCTION

### 1.1 Scope

The analytical method presented is used for the determination of parent milbemectin ( $M.A_3$  and  $M.A_4$ ), and its potential degradation products including the photoproducts ( $8,9Z-M.A_3$  and  $8,9Z-M.A_4$ ), hydroxy metabolites ( $27\text{-hydroxy-}M.A_3$  and  $27\text{-hydroxy-}M.A_4$ ) and keto metabolites ( $27\text{-keto-}M.A_3$  and  $27\text{-keto-}M.A_4$ )

The method may be used for the simultaneous determination of the combined residue of  $M.A_3$  and  $8,9Z-M.A_3$  and the combined residue of  $M.A_4$  and  $8,9Z-M.A_4$ . This is a direct result of the conversion of  $M.A_3$  and  $8,9Z-M.A_3$  to the same anhydro- $M.A_3$  compound resulting from the derivatization procedure employed to generate compounds that can be detected and quantitated. Likewise, the combined residue of  $M.A_4$  and  $8,9Z-M.A_4$  form the same anhydro- $M.A_4$  derivative via this derivatization procedure and, thus, can also be quantitated as a single combined residue.

Alternatively,  $M.A_3$ ,  $M.A_4$ ,  $8,9Z-M.A_3$  and  $8,9Z-M.A_4$  may be determined as individual components by incorporation of a fractionation step prior to derivatization of a composite extract. The resulting individual fractions would then be derivatized and analyzed separately for quantitation of the respective anhydro derivatives.

The procedure used for extraction of  $M.A_3$ ,  $M.A_4$ ,  $8,9Z-M.A_3$  and  $8,9Z-M.A_4$  from soil concurrently extracts the hydroxy metabolites ( $27\text{-hydroxy-}M.A_3$  and  $27\text{-hydroxy-}M.A_4$ ) and keto metabolites ( $27\text{-keto-}M.A_3$  and  $27\text{-keto-}M.A_4$ ). Methodology for the analysis of  $27\text{-keto-}M.A_3$  and  $27\text{-keto-}M.A_4$  consists of a direct chromatographic method since the keto metabolites do not form anhydro derivatives using the derivatization procedure presented in this method. Analysis of the  $27\text{-hydroxy-}M.A_3$  and  $27\text{-hydroxy-}M.A_4$  follow a similar methodology scheme to that of the parent compounds and their photoproducts, i.e., the same initial extraction, solid-phase clean up, and derivatization. The hydroxy anhydro derivatives are analyzed using chromatographic conditions used for either parent and/or photoproduct anhydro derivatives.

Limits of quantitation (LOQ) achievable using the method presented below are summarized as follows

- 1.00  $\mu\text{g a.i./Kg}$  for  $M.A_3$  and  $M.A_4$
- 1.12  $\mu\text{g a.i./Kg}$  for  $8,9Z-M.A_3$  and  $8,9Z-M.A_4$
- 1.00  $\mu\text{g a.i./Kg}$  for  $27\text{-hydroxy-}M.A_3$  and  $27\text{-hydroxy-}M.A_4$
- 2.00  $\mu\text{g a.i./Kg}$  for  $27\text{-keto-}M.A_3$  and  $27\text{-keto-}M.A_4$

### 1.2 Principle


Milbemectin, a mixture of  $M.A_3$  and  $M.A_4$ ; the photoproducts of milbemectin,  $8,9Z-M.A_3$  and  $8,9Z-M.A_4$ ; and the potential degradation products  $27\text{-hydroxy-}M.A_3$ ,  $27\text{-hydroxy-}M.A_4$ ,  $27\text{-keto-}M.A_3$  and  $27\text{-keto-}M.A_4$  are extracted from soil using methanol. The extracts are cleaned up using a  $C_{18}$  solid-phase extraction column. Residues are eluted from the column, the eluant evaporated to dryness and the residue reconstituted in

hexane/isopropanol. At this point, an aliquot of reconstituted eluant can be fractionated using high performance liquid chromatography (HPLC) to separate M.A<sub>3</sub>/M.A<sub>4</sub> from 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub>. Fractionation provides individual fractions, which can be dried and derivatized, thus allowing quantitation of the individual components.

For analysis of total<sup>1</sup> M.A<sub>3</sub> (including any 8,9Z-M.A<sub>3</sub> present), total<sup>1</sup> M.A<sub>4</sub> (including any 8,9Z-M.A<sub>4</sub> present), 27-hydroxy-M.A<sub>3</sub>, 27-hydroxy-M.A<sub>4</sub>, 27-keto-M.A<sub>3</sub>, and 27-keto-M.A<sub>4</sub>, a volumetric aliquot of the reconstituted eluant from the C<sub>18</sub> SPE column is transferred to a centrifuge tube and the residue subsequently derivatized to the fluorescence-active anhydro-M.A<sub>3</sub>, anhydro-M.A<sub>4</sub>, anhydro-27-hydroxy-M.A<sub>3</sub> and anhydro-27-hydroxy-M.A<sub>4</sub> derivatives. An aliquot of the remaining column eluant is evaporated to dryness, redissolved in hexane/isopropanol and analyzed directly for 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub> using gradient-elution, normal phase HPLC with ultraviolet absorption (UV) detection. As noted above, this procedure generates results for total M.A<sub>3</sub> that includes M.A<sub>3</sub> and its corresponding photoproduct, 8,9Z-M.A<sub>3</sub>, since they form the same anhydro-M.A<sub>3</sub> derivative. Likewise, the above procedure provides results for total M.A<sub>4</sub> that includes its corresponding photoproduct, 8,9Z-M.A<sub>4</sub>, since they form the same anhydro-M.A<sub>4</sub> derivative. Also, as noted above, the method may be altered to obtain separate values for parent and photoproduct components. To determine separate concentrations for these components, an aliquot of the solid-phase column eluant must be fractionated prior to performing the derivatization with collection of separate fractions for the parent and photoproduct moieties.

All anhydro derivatives, to include anhydro-M.A<sub>3</sub> (arising from either M.A<sub>3</sub> or 8,9Z-M.A<sub>3</sub>), anhydro-M.A<sub>4</sub> (arising from either M.A<sub>4</sub> or 8,9Z-M.A<sub>4</sub>), anhydro-27-hydroxy-M.A<sub>3</sub> and anhydro-27-hydroxy-M.A<sub>4</sub>, are analyzed chromatographically using gradient-elution, reversed-phase HPLC with fluorescence detection.

The derivatization reaction is illustrated below using milbemectin components M.A<sub>3</sub> and M.A<sub>4</sub> as examples. The analogous reaction occurs for the photoproducts 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>, and the hydroxy metabolites, 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>.

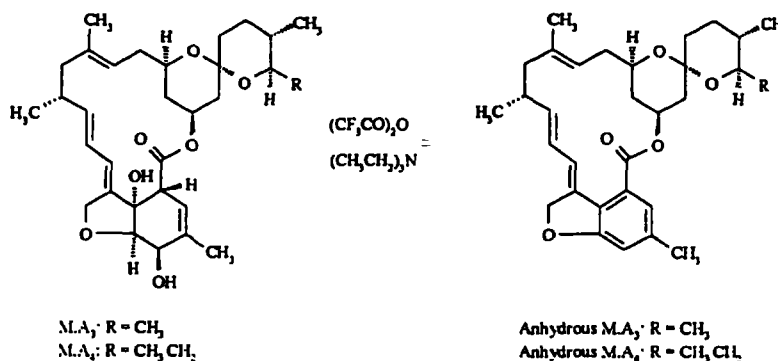

<sup>1</sup> The reference to "total M.A<sub>3</sub>" or "total M.A<sub>4</sub>" within this report refers to the combined residue of M.A<sub>3</sub> and its photoproduct 8,9Z-M.A<sub>3</sub> or M.A<sub>4</sub> and its photoproduct 8,9Z-M.A<sub>4</sub>, respectively.



WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 8 -

It is apparent from the above discussion that the methodology presented offers several options for analysis depending upon the moieties for which quantitative information is desired. To aid the analyst, the method is illustrated diagrammatically in Figure 3 with reference to the appropriate analytical sections.

## 2.0 APPARATUS

Note: Apparatus listed may be replaced with equivalent apparatus or other apparatus serving the same function if verification supports such substitutions.

- 2.1 Syringes, assorted sizes
- 2.2 Volumetric flasks, 100 and 200 mL, Class A, Kimax
- 2.3 Beakers, Kimax, 150 - 250 mL
- 2.4 Volumetric pipets, 10 and 20 mL, Class A, Kimax
- 2.5 Centrifuge tubes, 15 mL, graduated to 0.10 mL, Kimax
- 2.6 Vacuum manifold, Vac Master 10, Catalog No. 12H016, Jones Chromatography
- 2.7 Centrifuge bottles, Pyrex, 150 mL, Catalog No. 1265-150, Corning
- 2.8 Sonic Dismembrator, Model 300 operated at 60% power, Fisher
- 2.9 Centrifuge, Centra-8R, IEC
- 2.10 Filtering apparatus, filter funnel, 140 mL, 60-mm disk diameter, EC with 24/40 male adapter, WDLF-0002, At-Mar Glass Company

## 3.0 REAGENTS

Note: Reagents may be replaced with equivalent reagents from alternative vendors if verification supports such substitutions.

- 3.1 Methanol, HPLC grade, (203 nm UV cutoff), Burdick and Jackson,
- 3.2 Acetonitrile, HPLC grade (188 nm UV cutoff), Burdick, and Jackson
- 3.3 Water, Barnstead NANOpure<sup>®</sup>, ASTM type II, 17.8 M $\Omega$ -cm
- 3.4 2-Propanol, HPLC grade, Fisher
- 3.5 n-Hexane, HPLC grade, Burdick & Jackson
- 3.6 Benzene, HPLC grade, Sigma Aldrich

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 9 -

- 3.7 Triethylamine (TEA), analytical reagent grade, Aldrich
- 3.8 Trifluoroacetic anhydride (TFAA), analytical reagent grade, Aldrich
- 3.9 Phosphoric acid, reagent grade, Mallinckrodt
- 3.10 C<sub>18</sub>-reversed phase solid-phase extraction columns, 1000 mg, Bakerbond
- 3.11 Glass fiber filters, GF/A, Whatman
- 3.12 Celite 545, Catalog No. 3371-05, J.T. Baker

## 4.0 REFERENCE STANDARDS IDENTIFICATION

M.A<sub>3</sub>: (10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-21,24-dihydroxy-5',6',11,13,22-pentamethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

M.A<sub>4</sub>: (10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-6'-ethyl-21,24-dihydroxy-5',11,13,22-tetramethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

8,9Z-M A<sub>3</sub>: (10E,14E,16Z,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-21,24-dihydroxy-5',6',11,13,22-pentamethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

8,9Z-M A<sub>4</sub>: (10E,14E,16Z,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-6'-ethyl-21,24-dihydroxy-5',11,13,22-tetramethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

27-hydroxy- M.A<sub>3</sub>: 27-hydroxy-(10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-21,24-dihydroxy-5',6',11,13,22-pentamethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one

27-hydroxy- M.A<sub>4</sub>: 27-hydroxy-(10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-6'-ethyl-21,24-dihydroxy-5',11,13,22-tetramethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

27-keto- M.A<sub>3</sub>: 27-keto-(10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-21,24-dihydroxy-5',6',11,13,22-pentamethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

27-keto- M.A<sub>4</sub>: 27-keto-(10E,14E,16E,22Z)-(1R,4S,5'S,6R,6'R,8R,13R,20R,21R,24S)-6'-ethyl-21,24-dihydroxy-5',11,13,22-tetramethyl-3,7,19-trioxacyclo[15.6.1.1<sup>4,8</sup>.O<sup>20,24</sup>]-pentacos-10,14,16,22-tetraene-6-spiro-2'-tetrahydropyran-2-one.

Note: Store reference standards frozen

## 5.0 ANALYTICAL STANDARD SOLUTION PREPARATION

5.1 Stock Solution Preparation - M.A<sub>3</sub> and M.A<sub>4</sub>

Stock solutions of milbemectin are prepared by accurately weighing 100 mg (corrected for purity of the analytical standard) of each of the standards, M.A<sub>3</sub> and M.A<sub>4</sub>, onto weigh paper, transferring to separate clean 100-mL class A volumetric flasks and adjusting to final volume using methanol. These stock solutions contain 1.00 mg a.i./mL of M.A<sub>3</sub> and M.A<sub>4</sub>, respectively. The solutions are then volumetrically diluted 10 mL to 100 mL using methanol to achieve 100-mg a.i./L stock solutions. Serial dilutions of 10 mL to 100 mL are successively made in methanol to prepare stock solutions containing 10.0 and 1.00 mg a.i./L of M.A<sub>3</sub> and M.A<sub>4</sub>. Combined stocks, containing both M.A<sub>3</sub> and M.A<sub>4</sub>, may also be prepared. These solutions are used for the preparation of standards and for fortification of recovery samples.

5.2 Stock Solution Preparation - 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>

Stock solutions of the photoproducts of milbemectin are prepared by accurately weighing 100 mg (corrected for purity of the analytical standards) of each of the standards, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>, onto weigh paper, transferring to clean separate 100-mL class A volumetric flasks and adjusting to final volume using methanol. These stock solutions contain 1.00 mg a.i./mL of 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>, respectively. The solutions are then volumetrically diluted 10 mL to 100 mL using methanol to achieve 100 mg a.i./L stock solutions. Serial dilutions of 10 mL to 100 mL are successively made in methanol to prepare stock solutions containing 10.0 and 1.00 mg a.i./L of 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>. Combined stocks, containing both 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub> may also be prepared. These solutions are used for fortification of the recovery samples for the photoproduct method validation trials.

5.3 Stock Solution Preparation - 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>

Stock solutions of the 27-hydroxy metabolites of milbemectin are prepared by accurately weighing 25 mg (corrected for purity of the analytical standards) of each of the standards, 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>, onto weigh paper, transferring to separate clean 50-mL class A volumetric flasks and adjusting to final volume using acetonitrile. These stock solutions contain 0.500 mg a.i./mL of 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>, respectively. The solutions are then volumetrically diluted 20 mL to 100 mL using acetonitrile to achieve 100-mg a.i./L stock solutions. Serial dilutions of 10 mL to 100 mL are successively made in acetonitrile to prepare stock solutions containing 10.0 and 1.00 mg a.i./L of 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>. Combined stocks, containing both 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub> may also be prepared. These solutions are used for fortification of the recovery samples for the 27-hydroxy metabolite method validation trials.

5.4 Stock Solution Preparation - 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>

Stock solutions of the 27-keto metabolites of milbemectin are prepared by accurately weighing 25 mg (corrected for purity of the analytical standards) of each of the standards,

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 11 -

27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, onto weigh paper, transferring to separate clean 50-mL class A volumetric flasks and adjusting to final volume using acetonitrile. These stock solutions contain 0.500 mg a.i./mL of 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, respectively. The solutions are then volumetrically diluted 20 mL to 100 mL using acetonitrile to achieve 100-mg a.i./L stock solutions. Serial dilutions of 10 mL to 100 mL are successively made in acetonitrile to prepare stock solutions containing 10.0 and 1.00 mg a.i./L of 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>. Combined stocks, containing both 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub> may also be prepared. These solutions are used for fortification of the recovery samples for the 27-keto metabolite method validation trials.

#### 5.5 Combined Stock Solution - M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>

Stock solutions containing the four components, M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>, are prepared for the total residue and fractionation portions of the method. Ten milliliters each of the 100 mg a.i./L methanolic stock solutions of M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub> are volumetrically combined in a 100-mL class A volumetric flask. The final volume is adjusted to 100 mL with methanol, resulting in a solution containing 10.0 mg a.i./L of each of the four components. From this solution, 10.0 mL is volumetrically transferred to a 100-mL volumetric flask partially filled with methanol and then brought to volume with methanol. This solution contains 1.00 mg a.i./L of each component.

#### 5.6 Calibration Standard Solutions of Anhydro-M.A<sub>3</sub> and Anhydro-M.A<sub>4</sub>

Since M.A<sub>3</sub> and M.A<sub>4</sub>, and the respective 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub> photoproducts form the same anhydro-M.A<sub>3</sub> and anhydro-M.A<sub>4</sub> derivatives, M.A<sub>3</sub> and M.A<sub>4</sub> are derivatized and used to quantitate both parent milbemectin and its photoproducts. Calibration standards are prepared by addition of aliquots of the M.A<sub>3</sub> and M.A<sub>4</sub> stock solutions prepared in Sections 5.1 and 5.2 to clean, dry, 15-mL, centrifuge tubes using a gas-tight syringe. The stock solutions are subsequently evaporated to dryness and derivatized as described in the Section 6.5, Derivatization.

Stock Concentration (mg a.i./L) <sup>1</sup>	Fortification Volume (μL)	Final Volume (mL)	Standard Concentration (μg a.i./L) <sup>2</sup>
0.100	35.0	10.0	0.350
0.100	100	10.0	1.00
0.100	250	10.0	2.50
0.100	350	10.0	3.50
0.100	500	10.0	5.00

<sup>1</sup> The stock concentration contains 0.100 mg a.i./L of M.A<sub>3</sub> and 0.100 mg a.i./L of M.A<sub>4</sub>.

<sup>2</sup> Solutions contain M.A<sub>3</sub> and M.A<sub>4</sub>, each at the concentration presented.

#### 5.7 Calibration Standard Solutions of Anhydro-27-Hydroxy-M.A<sub>3</sub> and Anhydro-27-Hydroxy-M.A<sub>4</sub>

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 12 -

Since 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub> form distinctly different anhydro-derivatives from M.A<sub>3</sub> and M.A<sub>4</sub>, namely anhydro-27-hydroxy-M.A<sub>3</sub> and anhydro-27-hydroxy-M.A<sub>4</sub>, a separate set of standards may be derivatized and used to quantitate 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>. Alternatively, calibration standards may be prepared by fortifying centrifuge tubes with the four components, M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>, derivatizing and diluting the mixture as a singular combined standard. This approach is preferable since standards containing the four anhydro-derivatives are simultaneously chromatographed and the use of a combined standard minimizes the number of injections, thereby lessening the total run time. All combined standards are fortified, evaporated, derivatized and diluted using the methodology described in Section 5.6.

#### 5.8 Calibration Standard Solutions of 27-Keto-M.A<sub>3</sub> and 27-Keto-M.A<sub>4</sub>

Using the derivatization procedures previously described, the 27-keto-moieties of milbemectin, 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, do not form the corresponding anhydro-27-keto-M.A<sub>3</sub> and anhydro-27-keto-M.A<sub>4</sub>. Therefore, a separate set of calibration standards containing 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub> must be prepared as follows:

Stock Concentration (mg a.i./L) <sup>1</sup>	Fortification Volume (μL)	Final Volume (mL)	Standard Concentration (μg a.i./L) <sup>2</sup>
10.0	40.0	100	4.00
10.0	100	100	10.0
10.0	150	100	15.0
10.0	200	100	20.0
10.0	250	100	25.0

<sup>1</sup> The stock concentration contains 10.0 mg a.i./L of 27-keto-M.A<sub>3</sub> and 10.0 mg a.i./L of 27-keto-M.A<sub>4</sub>.

<sup>2</sup> Solutions contain 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, each at the concentration presented.

## 6.0 ANALYTICAL METHOD

### 6.1 Sample Preparation

- 6.1.1 Using an appropriate technique, mix soil in order to ensure a homogenous sample.
- 6.1.2 Weigh 25 grams of soil into a tared, pre-labeled, 150-mL, glass centrifuge bottle. Record the weight.

### 6.2 Laboratory Fortification

- 6.2.1 For the preparation of procedural recovery samples, soil is fortified as appropriate using stock solutions prepared in Sections 5.1 through 5.5.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

- 13-

- 6.2.2 For fractionation of M.A<sub>3</sub> and M.A<sub>4</sub> from 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>, fortify the soil samples with combined stock solutions of the respective components. Typical fortification levels for the fractionation method (parent/photoproduct) are 2.00, 10.0, 50.0 and 100 µg a.i./Kg
- 6.2.3 For analysis of M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub>, 27-hydroxy-M.A<sub>4</sub>, 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, fortify the soil samples with combined stock solutions of the six components. Typical fortification levels for the combined method are 2.00, 10.0 and 100 µg a.i./Kg for M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub>, 27-hydroxy-M.A<sub>4</sub> and 5.00, 10.0 and 100 µg a.i./Kg for 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>.

### 6.3 Sample Extraction

- 6.3.1 To the centrifuge bottle containing the sample, add 50 mL of methanol using a graduated cylinder or tiltapet.
- 6.3.2 Using a sonic disruptor set at 60% power, disrupt the sample for approximately five minutes. The sample at this stage should be converted to fine particles. If another sample needs to be immediately extracted, rinse the disruptor horn with a small quantity of methanol and collect the rinsate in the centrifuge bottle containing the sample.
- 6.3.3 Cap the centrifuge bottle and centrifuge the bottle at a setting of 2500 rpm for five minutes.
- 6.3.4 Prepare a filtering aid by adding 10 grams of Celite 545 to a glass beaker and adding sufficient reagent grade water to make the mixture flowable. Using vacuum, pour this mixture into a filtering apparatus with sintered glass frit (or Buchner funnel) fitted with a glass fiber filter. Ensure that the filtering aid forms an even, homogenous bed with a depth of approximately 5 to 10 mm of Celite.
- 6.3.5 Gently remove the sample from the centrifuge, ensuring that the solids are not disturbed. Vacuum filter the supernatant through the filtering aid and collect the filtrate in a labeled, 200-mL, glass, Class A, volumetric flask.
- 6.3.6 Add 50 mL of methanol to the centrifuge bottle and return to the sonic disruptor. Repeat steps 6.3.2, 6.3.3 and 6.3.5 combining the filtrate with the initial filtrate in the volumetric flask. The same filter aid prepared and used in step 6.3.4 for filtering the initial extract should also be used for the second extract.
- 6.3.7 Rinse the probe with 5 to 10 mL of methanol and collect the rinsate in the centrifuge bottle. Filter the rinsate through the Celite and collect in the same volumetric flask.

Note: The total volume of methanol in the flask must not exceed 135 mL or breakthrough during the subsequent solid-phase extraction step may occur.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

- 14 -

6.3.8 To the volumetric flask, add 65 mL of reagent grade water. Adjust the final volume to a final volume of 200-mL with methanol. Invert several times to ensure homogeneity. This solution will consist of 70% or less methanol and approximately 30% or greater water.

#### 6.4 Solid-Phase Extraction

6.4.1 For each sample, prepare a Bakerbond C<sub>18</sub> SPE (1000 mg) column by rinsing the column with 2 to 3 column volumes of methanol followed by 2 to 3 column volumes of reagent-grade water. Ensure that the column does not go to dryness until step 6.4.4.

6.4.2 Using a Class A volumetric pipet, transfer a portion of the extract (from step 6.4.1) to the column. In order to achieve the specified LOQs, a 20-mL aliquot of the extract is required.

6.4.3 Using a low vacuum, allow the sample to pass through the C<sub>18</sub> column. When the sample just reaches the top of the C<sub>18</sub> packing, add 5 mL of CH<sub>3</sub>OH/H<sub>2</sub>O, 70:30 (v:v) and continue elution. Successively add two additional 5 mL rinses of CH<sub>3</sub>OH/H<sub>2</sub>O, 70:30 (v:v). Discard all eluants.

6.4.4 Allow the column to elute to dryness and leave under vacuum for approximately 15 minutes in order to remove residual water.

Note: Increase vacuum as necessary. It is important that the column be allowed to dry completely under vacuum since residual water will slow the subsequent evaporation process and if not removed, will inhibit the formation of the anhydro-derivatives.

6.4.5 Volumetrically add 5 mL of 100% methanol to the column and elute into a 15-mL, graduated, centrifuge tube.

6.4.6 Place the centrifuge tube in a waterbath maintained at approximately 40°C and evaporate the extract to dryness under a gentle stream of nitrogen.

Note. Ensure that no water is present in the centrifuge tube prior to derivatization. If a small quantity of water remains following evaporation, it may be azeotroped by addition of absolute ethanol followed by evaporation under nitrogen.

#### 6.5 Separation of Milbemectin and its Potential Degradation Products

6.5.1 For analysis of total M<sub>A3</sub> and M<sub>A4</sub> residues, proceed to section 6.6.

6.5.2 For fractionation of M<sub>A3</sub>/M<sub>A4</sub> from 8.9Z-M<sub>A3</sub>/8.9Z-M<sub>A4</sub>, proceed to section 6.7

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 15-

6.5.2 For analysis of M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub>, 27-hydroxy-M.A<sub>4</sub>, 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, proceed to section 6.8.

#### 6.6 Derivatization

6.6.1 To the centrifuge tube, add 1.00 mL of 0.5M triethylamine (TEA) in benzene and 100 µL of trifluoroacetic anhydride. Vortex the tube for a few seconds to ensure mixing. Place the uncapped tube in a water bath maintained at approximately 40°C and allow the reaction to proceed for thirty minutes.

6.6.2 Remove the tube from the water bath and add 50 µL of 100% triethylamine. Vortex the solution for a few seconds.

6.6.3 Return the tube to the 40°C water bath and evaporate the solution under a gentle stream of nitrogen.

Note: The solution will not evaporate to dryness due to the presence of triethylamine. The final volume should be approximately 0.2 mL.

6.6.4 Adjust the final volume of each sample (using the graduations on the centrifuge tube) with CH<sub>3</sub>CN:H<sub>2</sub>O, 95:5 (v:v). Vortex the solution for a few seconds and transfer a portion to an autosampler vial for analysis.

Note: The final volume will be dictated by the residue level present. In order to achieve an LOQ of 1.00 µg a.i./Kg, a final volume of 5.00 mL is required when the low standard is 0.500µg a.i./L.

#### 6.7 Separation of M.A<sub>3</sub>/M.A<sub>4</sub> and 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub>

An aliquot of the reconstituted residue resulting following drying of the eluent of the solid-phase column clean up may be fractionated to separate parent milbemectin and its photoproducts using the following procedure.

6.7.1 Dissolve the residue (from step 6.4.6) in a centrifuge tube by volumetrically adding 1.00 mL of a solution of n-hexane:2-propanol, 99:1, (v:v). Vortex the centrifuge tube for approximately 30 seconds.

6.7.2 If determination of the total milbemectin residue is desired, an aliquot of the sample may be taken for subsequent derivatization/analysis with an additional aliquot retained for determination of the individual components. For analysis of total residues, volumetrically transfer 250 µL of the solution to a 15-mL graduated centrifuge tube, evaporate to dryness under a gentle stream of nitrogen and derivatize by following steps 6.6.1 through 6.6.4.

6.7.3 To separate M.A<sub>3</sub>/M.A<sub>4</sub> from the photoproducts 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub>, transfer an aliquot of the reconstituted dried eluent to an HPLC autosampler vial for injection and fractionation using normal phase HPLC.



WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

- 16-

- 6.7.4 Ensure that the HPLC column is equilibrated and that the instrument is fully operational. Chromatographic conditions for the separation of M.A<sub>3</sub>/M.A<sub>4</sub> from 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub> are presented in Section 7 of this method.
- 6.7.5 Prepare a reference standard in n-hexane:2-propanol, 99:1, (v:v) that contains M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub> at such a level that the components can be detected by UV using the supplied operating conditions. Inject this standard to ensure retention time stability.
- 6.7.6 Identify the retention times of the four components. Using the chromatographic conditions summarized in Section 7 of this method, the retention time of the photoproducts should be in the range of 8 to 12 minutes; the parent components in the range of 12 to 16 minutes. Retention ranges must be confirmed daily.
- 6.7.7 Inject samples and collect fractions as follows. Define the time range for the photoproducts and parent components such that the column eluant is collected well before and after the peaks of interest. Samples may be collected manually or with the use of an automated fraction collector.
- 6.7.8 Inject the reference standard after chromatographing the sample set to ensure retention time stability was maintained.
- 6.7.9 Evaporate the eluant fractions to dryness under a gentle stream of nitrogen. Proceed with derivatization, as presented in Sections 6.6.1 through 6.6.4.
- 6.8 Fractionation of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub>/27-hydroxy-M.A<sub>4</sub> and 27-keto-M.A<sub>3</sub>/27-keto-M.A<sub>4</sub>

An aliquot of the reconstituted residue resulting following drying of the eluant of the solid-phase column clean up may be fractionated to separate parent milbemectin, the 27-hydroxy metabolites and the 27-keto metabolites using the following procedure.

- 6.8.1 Dissolve the residue (from step 6.4.6) in the centrifuge tube by volumetrically adding 5.00 mL of a solution of n-hexane:2-propanol, 99:1, (v:v). Vortex the centrifuge tube for approximately 30 seconds.
- 6.8.2 For analysis of M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub>, volumetrically transfer 1000 µL of the reconstituted solution to a second graduated centrifuge tube.
- 6.8.3 Evaporate the 1000 µL of solution to dryness under a gentle stream of nitrogen. Proceed with derivatization, as presented in Sections 6.6.1 through 6.6.4.
- 6.8.4 For analysis of 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub>, evaporate the remaining 4 mL of extract to dryness under a gentle stream of nitrogen. Reconstitute to volume with n-hexane:2-propanol, 99:1 (v:v) and submit samples to analysis

WILDLIFE INTERNATIONAL LTD.

PROJECT NO. 420C-106

- 17 -

by normal-phase HPLC. Chromatographic conditions are presented in section 7.3.

### 6.9 HPLC Determination

Concentrations of residues that have been converted via derivatization to their anhydro derivatives are determined using HPLC with fluorescence detection. Concentrations of the keto metabolites are determined using HPLC with UV detection. Complete chromatographic conditions are presented in Section 7.0.

## 7.0 CHROMATOGRAPHIC CONDITIONS

### 7.1 HPLC Conditions for the analysis of Anhydro-M.A<sub>3</sub>, Anhydro-M.A<sub>4</sub>, Hydroxy M.A<sub>3</sub> and Hydroxy M.A<sub>4</sub> Derivatives

Instrument: Hewlett Packard Model 1090 High Performance Liquid Chromatograph  
 Detector: Jasco Model FP-920 Fluorescence Detector  
 Analytical Column: YMC-Pack ODS-AM 250 mm x 4.6 mm, 3- $\mu$ m particle size  
 Flow Rate: 1.0 mL/min  
 Column Temperature: 40°C  
 Mobile Phase Composition:

Solvent A:  
20% CH<sub>3</sub>CN: 80% H<sub>2</sub>O: 0.1% H<sub>3</sub>PO<sub>4</sub>

Solvent B:  
98%CH<sub>3</sub>CN: 2% H<sub>2</sub>O: 0.1% H<sub>3</sub>PO<sub>4</sub>

Time (min)	%A	%B
0.01	45	55
1.00	45	55
15.00	0	100
25.00	0	100
25.10	45	55
30.00	45	55

Injection Volume: 100  $\mu$ L  
 Retention Times of Anhydro-M.A<sub>3</sub> Approximately 20.8 minutes  
 Anhydro M.A<sub>4</sub> Approximately 22.5 minutes  
 27-Hydroxy M.A<sub>3</sub> Approximately 15.5 minutes  
 27-Hydroxy M.A<sub>4</sub> Approximately 16.8 minutes  
 Excitation/Emission Wavelength: 360/460 nm

### 7.2 Normal Phase HPLC Conditions for the Fractionation of 8,9Z-M.A<sub>3</sub>, 8,9Z-M.A<sub>4</sub>, M.A<sub>3</sub> and M.A<sub>4</sub>

Instrument: Hewlett-Packard Model 1090 High Performance Liquid Chromatograph  
 Detector: Diode Array Detector  
 Analytical Column: Phenomenex LUNA Silica 250 mm x 4.6 mm, 5- $\mu$ m particle size

WILDLIFE INTERNATIONAL LTD.

PROJECT NO. 420C-106

- 18 -

Flow Rate: 2.0 mL/min  
 Column Temperature: 40°C  
 Mobile Phase Composition: 99:1, Hexane : 2-Propanol (v:v)  
 Injection Volume: 250 µL  
 Retention Times of M.A<sub>3</sub>: Approximately 20.7 minutes  
 M.A<sub>4</sub>: Approximately 18.3 minutes  
 8,9Z-M.A<sub>3</sub>: Approximately 13.6 minutes  
 8,9Z-M.A<sub>4</sub>: Approximately 13.0 minutes  
 Primary Analytical Wavelength: 244 nm

### 7.3 Isocratic HPLC Operational Parameters For The Analysis of Anhydro-M.A<sub>3</sub> and Anhydro-M.A<sub>4</sub>

Instrument: Hewlett Packard Model 1090 High Performance Liquid Chromatograph  
 Detector: Jasco Model FP-920 Fluorescence Detector  
 Analytical Column: HP ODS Hypersil 200 x 4.6 mm, 5 µm  
 Flow Rate: 1.0 mL/min  
 Column Temperature: 40°C  
 Mobile Phase Composition: 95:5 Acetonitrile : Water (v:v)  
 Injection Volume: 100 µL  
 Retention Times of Anhydro-M.A<sub>3</sub>: Approximately 7.4 minutes  
 Anhydro-M.A<sub>4</sub>: Approximately 8.6 minutes  
 Excitation/Emission Wavelength: 360/460 nm

### 7.4 Normal Phase HPLC Conditions for the Analysis of 27Keto-M.A<sub>3</sub> and 27Keto-M.A<sub>4</sub>

Instrument: Waters Model 1090 High Performance Liquid Chromatograph  
 Detector: Waters Model 486 Detector  
 Analytical Column: YMC-Silica 250 mm x 2.0 mm, 5-µm particle size  
 Flow Rate: 1.0 mL/min  
 Column Temperature: 40°C  
 Mobile Phase Composition: Solvent A: 99% Hexane 1% IPA  
 Solvent B: 95% Hexane 5% IPA

Time (min)	%A	%B
0.01	100	0
1.00	100	0
14.00	0	100
14.10	100	0
18.00	100	0

Injection Volume: 200 µL  
 Retention Times of 27Keto-M.A<sub>3</sub>: Approximately 10.0 minutes  
 27Keto-M.A<sub>4</sub>: Approximately 9.5 minutes  
 Primary Analytical Wavelength: 280 nm

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 19 -

## 8.0 CALIBRATION

Calibration standards containing M.A<sub>3</sub>, M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub> and 27-hydroxy-M.A<sub>4</sub> were prepared by fortifying dry 15 mL graduated centrifuge, evaporating the solvent to dryness under a gentle stream of nitrogen and derivatizing according to sections 6.6.1-6.6.4. These standards yielded concentrations ranging from 0.350-5.00 µg a.i./L. Linear regression equations were generated using the respective concentrations of the calibration standards versus peak area responses of the anhydro-M.A<sub>3</sub>, anhydro-M.A<sub>4</sub>, anhydro-27-hydroxy-M.A<sub>3</sub>, and anhydro-27-hydroxy-M.A<sub>4</sub> derivatives, as follows

$$\text{Concentration } (\mu\text{g a.i./Kg}) = \frac{(\text{Peak area} - \text{y-intercept}) \times \text{dilution factor}}{\text{Slope}}$$

Four representative calibration curves for the parent and the 27-hydroxy components for a standard range of 0.350-5.00 µg a.i./L are presented in Figure 4, 5, 6 and 7

Calibration standards containing 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub> at concentrations ranging from 4.00-25.0 µg a.i./L were fortified directly into 100 mL volumetric flasks containing 99.1 hexane:IPA (v/v). Two representative calibration curves for the 27-keto derivatives are presented in figures 8 and 9

## 9.0 CALCULATIONS

The concentration of M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub>, 8,9Z-M.A<sub>4</sub>, 27-hydroxy-M.A<sub>3</sub>, 27-hydroxy-M.A<sub>4</sub>, 27-keto-M.A<sub>3</sub> and 27-keto-M.A<sub>4</sub> (milbemectin component) in samples is determined by substitution of the peak area response into the regression equation obtained from calibration standards

$$\text{Milbemectin component } (\mu\text{g a.i./Kg}) = [(\text{Peak Area} - \text{Y-Intercept})/\text{Slope}] \times \text{Dilution Factor}$$

where the Dilution Factor adjusts for the dilution and/or concentration of sample by the method

$$\text{Dilution Factor} = \frac{V_{12}}{V_{11}} \times \frac{V_{11}}{M_1}$$

Where:

M<sub>1</sub> = Initial mass of soil extracted

V<sub>11</sub> = Final extracted volume

V<sub>12</sub> = Volume of extract on the c<sub>18</sub> column

V<sub>12</sub> = Final volume

Recoveries for fortified samples are calculated as follows.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.. 420C-106

- 20 -

$$\text{Percent Recovery} = \frac{\text{Measured milbemectin component } (\mu\text{g a.i./Kg})}{\text{Nominal milbemectin component } (\mu\text{g a.i./Kg})} \times 100$$

An example calculation using a sample fortified at 10.0  $\mu\text{g a.i./Kg}$  with M.A<sub>3</sub> follows:

where      Peak Area = 517 44684  
               Y-Intercept = 7.44679  
               Slope = 221.6711  
               Dilution factor = 4.00

$$\text{Dilution Factor} = \frac{5.00}{10.0} \times \frac{200}{25.0} = 4.00$$

$$\begin{aligned} \text{Milbemectin in sample } (\mu\text{g a.i./Kg}) &= \frac{(517.44684 - 7.44679)}{221.6711} \times 4.00 \\ &= 9.20 \mu\text{g a.i./Kg} \end{aligned}$$

$$\text{Percent Recovery} = \frac{9.20 \mu\text{g a.i./Kg milbemectin}}{10.0 \mu\text{g a.i./Kg milbemectin}} \times 100$$

Percent Recovery = 92.0%

Representative chromatograms of low and high-level calibration standards of anhydro-M.A<sub>3</sub> and anhydro-M.A<sub>4</sub> are shown in figures 12 and 13 respectively.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 21 -

Table 1

Typical Gradient HPLC Conditions For The Analysis of M.A<sub>3</sub>, M.A<sub>4</sub>,  
27-Hydroxy M.A<sub>3</sub> and 27-Hydroxy M.A<sub>4</sub>

**INSTRUMENT:** Waters Model 2690 High Performance Liquid Chromatograph with a Jasco FP-920 Detector

**DATA ACQUISITION:** Millennium version 2.15

**ANALYTICAL COLUMN:** YMC-Pack ODS-AM Column (250 mm x 4.6 mm, 3 µm particle size)

**OVEN TEMPERATURE:** 40°C

**FLOW RATE:** 1.0 mL/minute

**MOBILE PHASE:** Solvent A: 20% Acetonitrile · 80% Water · 0.1% Phosphoric Acid  
Solvent B: 98% Acetonitrile · 2% Water · 0.1% Phosphoric Acid

<u>Time (min)</u>	<u>%A</u>	<u>%B</u>
0.01	45	55
1.00	45	55
15.00	0	100
25.00	0	100
25.10	45	55
30.00	45	55

**INJECTION VOLUME:** 100 µL

**RETENTION TIME:** M.A<sub>3</sub> : ~ 20.8  
M.A<sub>4</sub> : ~ 22.5  
27-Hydroxy M.A<sub>3</sub> : ~ 15.5  
27-Hydroxy M.A<sub>4</sub> : ~ 16.8

**FLUORESCENCE WAVELENGTH:** Excitation λ: 360 nm  
Emission λ: 460 nm

WILDLIFE INTERNATIONAL LTD.

PROJECT NO . 420C-106

- 22 -

Table 2

Typical Normal Phase HPLC Conditions For The Fractionation of 8,9Z-M.A<sub>3</sub>, 8,9Z-M.A<sub>4</sub>,  
M A<sub>3</sub> and M.A<sub>4</sub>

INSTRUMENT:	Hewlett-Packard Model 1090 High Performance Liquid Chromatograph with a Diode Array Detector
DATA ACQUISITION	HP Chemstation version 8.04
ANALYTICAL COLUMN:	Phenomenex LUNA Silica (250 mm x 4.6 mm, 5 µm particle size)
OVEN TEMPERATURE.	40°C
FLOW RATE.	2.0 mL/minute
MOBILE PHASE:	99.1, Hexane :2-Propanol (v:v)
INJECTION VOLUME:	250 µL
RETENTION TIME:	M.A <sub>3</sub> : ~20.7 minutes 8,9Z-M.A <sub>3</sub> : ~13.6 minutes M.A <sub>4</sub> ~ 18.3 minutes 8,9Z-M.A <sub>4</sub> : ~13 minutes
PRIMARY ANALYTICAL WAVELENGTH:	244nm

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 23 -

Table 3

Typical Normal Phase HPLC Conditions For The Analysis of 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

INSTRUMENT Waters Model 1090 High Performance Liquid Chromatograph  
with a Waters Model 486 Detector

DATA ACQUISITION. Millennium version 2.15

ANALYTICAL COLUMN YMC Silica (250 mm x 2.0 mm, 5 µm particle size)

OVEN TEMPERATURE. 40°C

FLOW RATE 1.0 mL/minute

MOBILE PHASE: Solvent A: 99% Hexane : 1% IPA  
Solvent B: 95% Hexane : 5% IPA

<u>Time (min)</u>	<u>%A</u>	<u>%B</u>
0.01	100	0
1.00	100	0
14.00	0	100
14.10	100	0
18.00	0	0

INJECTION VOLUME: 200 µL

RETENTION TIME: 27-Keto M.A<sub>3</sub> : ~10.0  
27-Keto M.A<sub>4</sub> : ~9.5

PRIMARY ANALYTICAL WAVELENGTH: 280 nm



WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 24 -

Table 4

Method Validation Recoveries of Fractionated M A<sub>3</sub> from New York Sandy Loam  
Fortified with M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-13	Reagent Blank	0.00	<1.12	--
VMAB-37	Matrix Blank	0.00	<1.12	--
VMAB-38	Matrix Blank	0.00	<1.12	--
VMAB-39	Matrix Blank	0.00	<1.12	--
VMAS-115	Matrix Fortification	2.00	1.49	74.3
VMAS-116	Matrix Fortification	2.00	1.62	80.8
VMAS-117	Matrix Fortification	2.00	1.72	86.1
VMAS-118	Matrix Fortification	10.0	8.80	88.0
VMAS-119	Matrix Fortification	10.0	8.82	88.2
VMAS-120	Matrix Fortification	10.0	8.88	88.8
VMAS-121	Matrix Fortification	50.0	44.0	88.0
VMAS-122	Matrix Fortification	50.0	41.3	82.6
VMAS-123	Matrix Fortification	50.0	41.5	83.0
VMAS-124	Matrix Fortification	100	77.7	77.7
VMAS-125	Matrix Fortification	100	81.0	81.0
VMAS-126	Matrix Fortification	100	79.4	79.4
Totals				Mean = 83.2
				S.D. = 4.71
				CV = 5.66
				N = 12

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.12 µg a.i./Kg as calculated from the product of the low standard (0.350 µg a.i./Kg) and the dilution factor of the matrix blanks (3.20).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.. 420C-106

- 25-

Table 5

Method Validation Recoveries of Fractionated M.A.<sub>4</sub> from New York Sandy Loam  
Fortified with M.A.<sub>3</sub>, M.A.<sub>4</sub>, 8,9Z-M.A.<sub>3</sub> and 8,9Z-M.A.<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-13	Reagent Blank	0.00	<1.12	--
VMAB-37	Matrix Blank	0.00	<1.12	--
VMAB-38	Matrix Blank	0.00	<1.12	--
VMAB-39	Matrix Blank	0.00	<1.12	--
VMAS-115	Matrix Fortification	2.00	1.50	75.2
VMAS-116	Matrix Fortification	2.00	1.68	83.8
VMAS-117	Matrix Fortification	2.00	1.77	88.7
VMAS-118	Matrix Fortification	10.0	9.07	90.7
VMAS-119	Matrix Fortification	10.0	9.09	90.9
VMAS-120	Matrix Fortification	10.0	9.12	91.2
VMAS-121	Matrix Fortification	50.0	45.3	90.6
VMAS-122	Matrix Fortification	50.0	42.6	85.2
VMAS-123	Matrix Fortification	50.0	43.0	86.0
VMAS-124	Matrix Fortification	100	81.0	81.0
VMAS-125	Matrix Fortification	100	83.7	83.7
VMAS-126	Matrix Fortification	100	82.8	82.8
Totals:				Mean = 85.8
				S.D. = 4.90
				CV = 5.71
				N = 12

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.12 µg a.i./Kg as calculated from the product of the low standard (0.350 µg a.i./Kg) and the dilution factor of the matrix blanks (3.20)

NOTE. Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly

WILDLIFE INTERNATIONAL LTD.

PROJECT NO. 420C-106

- 26 -

Table 6

Method Validation Recoveries of Fractionated 8,9Z-M.A<sub>3</sub> from New York Sandy Loam  
Fortified with M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-13	Reagent Blank	0.00	<1.12	--
VMAB-37	Matrix Blank	0.00	<1.12	--
VMAB-38	Matrix Blank	0.00	<1.12	--
VMAB-39	Matrix Blank	0.00	<1.12	--
VMAS-115	Matrix Fortification	2.00	2.43	122
VMAS-116	Matrix Fortification	2.00	1.91	95.6
VMAS-117	Matrix Fortification	2.00	2.00	100
VMAS-118	Matrix Fortification	10.0	10.23	102
VMAS-119	Matrix Fortification	10.0	9.25	92.5
VMAS-120	Matrix Fortification	10.0	9.26	92.6
VMAS-121	Matrix Fortification	50.0	44.6	89.2
VMAS-122	Matrix Fortification	50.0	44.7	89.3
VMAS-123	Matrix Fortification	50.0	47.2	94.4
VMAS-124	Matrix Fortification	100	84.0	84.0
VMAS-125	Matrix Fortification	100	88.2	88.2
VMAS-126	Matrix Fortification	100	88.4	88.4
Totals:				Mean = 94.9
				S.D. = 9.91
				CV = 10.44
				N = 12

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.12 µg a.i./Kg as calculated from the product of the low standard (0.350 µg a.i./Kg) and the dilution factor of the matrix blanks (3.20).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 27 -

Table 7

Method Validation Recoveries of Fractionated 8,9Z-M.A<sub>4</sub> from New York Sandy Loam  
Fortified with M A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M A<sub>3</sub> and 8,9Z-M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-13	Reagent Blank	0.00	<1.12	--
VMAB-37	Matrix Blank	0.00	<1.12	--
VMAB-38	Matrix Blank	0.00	<1.12	--
VMAB-39	Matrix Blank	0.00	<1.12	--
VMAS-115	Matrix Fortification	2.00	2.02	101
VMAS-116	Matrix Fortification	2.00	1.82	91.0
VMAS-117	Matrix Fortification	2.00	1.88	93.8
VMAS-118	Matrix Fortification	10.0	9.67	96.7
VMAS-119	Matrix Fortification	10.0	8.77	87.7
VMAS-120	Matrix Fortification	10.0	8.78	87.8
VMAS-121	Matrix Fortification	50.0	42.2	84.3
VMAS-122	Matrix Fortification	50.0	42.2	84.4
VMAS-123	Matrix Fortification	50.0	44.7	89.3
VMAS-124	Matrix Fortification	100	79.6	79.6
VMAS-125	Matrix Fortification	100	83.4	83.4
VMAS-126	Matrix Fortification	100	83.4	83.4
Totals:				Mean = 88.5
				S.D. = 6.20
				CV = 7.01
				N = 12

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.12 µg a.i./Kg as calculated from the product of the low standard (0.350 µg a.i./Kg) and the dilution factor of the matrix blanks (3.20).

NOTE. Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly

WILDLIFE INTERNATIONAL LTD

PROJECT NO.. 420C-106

- 28 -

Table 8

Method Validation Recoveries of M.A<sub>3</sub> from New York Sandy Loam  
 Fortified with M A<sub>3</sub>, M A<sub>4</sub>, 27-Hydroxy M.A<sub>3</sub>, 27-Hydroxy M.A<sub>4</sub>, 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<1.00	--
VMAB-88	Matrix Blank	0.00	<1.00	--
VMAB-89	Matrix Blank	0.00	<1.00	--
VMAB-90	Matrix Blank	0.00	<1.00	--
VMAS-280	Matrix Fortification	2.00	1.59	79.7
VMAS-281	Matrix Fortification	2.00	1.89	94.3
VMAS-282	Matrix Fortification	2.00	1.50	75.1
VMAS-283	Matrix Fortification	10.0	8.15	81.5
VMAS-284	Matrix Fortification	10.0	9.15	91.5
VMAS-285	Matrix Fortification	10.0	9.37	93.7
VMAS-286	Matrix Fortification	100	92.2	92.2
VMAS-287	Matrix Fortification	100	83.4	83.4
VMAS-288	Matrix Fortification	100	90.2	90.2
			Totals:	Mean = 86.9
				S.D = 7.01
				CV = 8.07
				N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.00 µg a.i./Kg as calculated from the product of the low standard (0.500 µg a.i./Kg) and the dilution factor of the matrix blanks (2.00).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 29 -

Table 9

Method Validation Recoveries of M.A<sub>4</sub> from New York Sandy Loam  
Fortified with M A<sub>3</sub>, M A<sub>4</sub>, 27-Hydroxy M.A<sub>3</sub>, 27-Hydroxy M.A<sub>4</sub>, 27-Keto M A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<1.00	--
VMAB-88	Matrix Blank	0.00	<1.00	--
VMAB-89	Matrix Blank	0.00	<1.00	--
VMAB-90	Matrix Blank	0.00	<1.00	--
VMAS-280	Matrix Fortification	2.00	1.63	81.4
VMAS-281	Matrix Fortification	2.00	1.91	95.4
VMAS-282	Matrix Fortification	2.00	1.54	76.9
VMAS-283	Matrix Fortification	10.0	8.31	83.1
VMAS-284	Matrix Fortification	10.0	9.26	92.6
VMAS-285	Matrix Fortification	10.0	9.46	94.6
VMAS-286	Matrix Fortification	100	92.5	92.5
VMAS-287	Matrix Fortification	100	83.9	83.9
VMAS-288	Matrix Fortification	100	91.1	91.1
			Totals.	Mean = 87.9
				S.D. = 6.66
				CV = 7.58
				N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.00 µg a.i./Kg as calculated from the product of the low standard (0.500 µg a.i./Kg) and the dilution factor of the matrix blanks (2.00).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO: 420C-106

- 30 -

Table 10

Method Validation Recoveries of 27-Hydroxy M.A<sub>3</sub> from New York Sandy Loam  
Fortified with M.A<sub>3</sub>, M.A<sub>4</sub>, 27-Hydroxy M.A<sub>3</sub>, 27-Hydroxy M.A<sub>4</sub>, 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<1.00	--
VMAB-88	Matrix Blank	0.00	<1.00	--
VMAB-89	Matrix Blank	0.00	<1.00	--
VMAB-90	Matrix Blank	0.00	<1.00	--
VMAS-280	Matrix Fortification	2.00	1.61	80.4
VMAS-281	Matrix Fortification	2.00	2.00	100
VMAS-282	Matrix Fortification	2.00	1.48	74.1
VMAS-283	Matrix Fortification	10.0	6.77	67.7
VMAS-284	Matrix Fortification	10.0	6.33	63.3
VMAS-285	Matrix Fortification	10.0	8.00	80.0
VMAS-286	Matrix Fortification	100	84.5	84.5
VMAS-287	Matrix Fortification	100	74.2	74.2
VMAS-288	Matrix Fortification	100	72.3	72.3
Totals:				Mean = 77.4
				S.D. = 10.7
				CV =
				13.82
				N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.00 µg a.i./Kg as calculated from the product of the low standard (0.500 µg a.i./Kg) and the dilution factor of the matrix blanks (2.00).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

- 31 -

Table 11

Method Validation Recoveries of 27-Hydroxy M A<sub>4</sub> from New York Sandy Loam  
Fortified with M.A<sub>3</sub>, M.A<sub>4</sub>, 27-Hydroxy M.A<sub>3</sub>, 27-Hydroxy M A<sub>4</sub>, 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<1.00	--
VMAB-88	Matrix Blank	0.00	<1.00	--
VMAB-89	Matrix Blank	0.00	<1.00	--
VMAB-90	Matrix Blank	0.00	<1.00	--
VMAS-280	Matrix Fortification	2.00	1.70	85.0
VMAS-281	Matrix Fortification	2.00	2.06	103
VMAS-282	Matrix Fortification	2.00	1.65	82.5
VMAS-283	Matrix Fortification	10.0	7.07	70.7
VMAS-284	Matrix Fortification	10.0	6.65	66.5
VMAS-285	Matrix Fortification	10.0	8.32	83.2
VMAS-286	Matrix Fortification	100	84.3	84.3
VMAS-287	Matrix Fortification	100	74.6	74.6
VMAS-288	Matrix Fortification	100	73.8	73.8
Totals				Mean = 80.4
				S.D. = 10.7
				CV = 13.31
				N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 1.00 µg a.i./Kg as calculated from the product of the low standard (0.500 µg a.i./Kg) and the dilution factor of the matrix blanks (2.00).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.



WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 32 -

Table 12

Method Validation Recoveries of 27-Keto M.A<sub>3</sub> from New York Sandy Loam  
 Fortified with M.A<sub>3</sub>, M.A<sub>4</sub>, 27-Hydroxy M.A<sub>3</sub>, 27-Hydroxy M.A<sub>4</sub>, 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<2.00	--
VMAB-88	Matrix Blank	0.00	<2.00	--
VMAB-89	Matrix Blank	0.00	<2.00	--
VMAB-90	Matrix Blank	0.00	<2.00	--
VMAS-280	Matrix Fortification	5.00	4.65	93.0
VMAS-281	Matrix Fortification	5.00	6.48	130
VMAS-282	Matrix Fortification	5.00	6.11	122
VMAS-283	Matrix Fortification	10.0	9.29	92.9
VMAS-284	Matrix Fortification	10.0	10.1	101
VMAS-285	Matrix Fortification	10.0	9.85	98.5
VMAS-286	Matrix Fortification	100	84.4	84.4
VMAS-287	Matrix Fortification	100	88.3	88.3
VMAS-288	Matrix Fortification	100	82.3	82.3
Totals:				Mean = 99.1
				S.D. = 16.4
				CV = 16.57
				N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 2.00 µg a.i./Kg as calculated from the product of the low standard (4.00 µg a.i./Kg) and the dilution factor of the matrix blanks (0.500).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

- 33 -

Table 13

Method Validation Recoveries of 27-Keto M A<sub>4</sub> from New York Sandy Loam  
 Fortified with M A<sub>3</sub>, M A<sub>4</sub>, 27-Hydroxy M A<sub>3</sub>, 27-Hydroxy M A<sub>4</sub>, 27-Keto M A<sub>3</sub> and 27-Keto M A<sub>4</sub>

Number (420C-106-)	Sample Type	Concentration (µg a.i./Kg)		Percent Recovery
		Fortified	Measured <sup>1</sup>	
VREB-30	Reagent Blank	0.00	<2.00	--
VMAB-88	Matrix Blank	0.00	<2.00	--
VMAB-89	Matrix Blank	0.00	<2.00	--
VMAB-90	Matrix Blank	0.00	<2.00	--
VMAS-280	Matrix Fortification	5.00	5.76	115
VMAS-281	Matrix Fortification	5.00	8.00	160
VMAS-282	Matrix Fortification	5.00	6.29	126
VMAS-283	Matrix Fortification	10.0	10.7	107
VMAS-284	Matrix Fortification	10.0	10.6	106
VMAS-285	Matrix Fortification	10.0	12.3	123
VMAS-286	Matrix Fortification	100	89.7	89.7
VMAS-287	Matrix Fortification	100	88.1	88.1
VMAS-288	Matrix Fortification	100	83.8	83.8
			Totals.	Mean = 99 S.D. = 16.4 CV = 16.57 N = 9

<sup>1</sup> The limit of quantitation (LOQ) was set at 2.00 µg a.i./Kg as calculated from the product of the low standard (4.00 µg a.i./Kg) and the dilution factor of the matrix blanks (0.500).

NOTE: Results were generated using Excel 5.0 in the full precision mode. Calculated values may differ slightly.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 34-

Figure 1 Analytical method flow chart for the analysis of Milbemectin and its metabolites in soil.

**METHOD OUTLINE FOR THE PROCESSING OF  
MILBEMECTIN AND ITS METABOLITES IN FIELD SOILS**

1. Weigh 25 grams of soil into a tared container or directly into a tared glass centrifuge tube.
2. If necessary, transfer the sample into a centrifuge tube. Fortify QC samples with the appropriate M.A<sub>3</sub>/M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub>/8,9-ZM.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> or 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> stock solutions.
3. To each centrifuge tube, add 50 mL of methanol.
4. Using a sonic disrupter, disrupt for approximately 5 minutes.
5. Centrifuge the mixture at 2500 rpm for 5 minutes.
6. Gently remove the sample from the centrifuge and vacuum filter the supernatant, using a glass fiber filter (containing approximately 10 grams of Celite 545 added as a slurry in water) placed on the sintered glass frit. Collect the filtrate in a 200-mL, class A volumetric flask.
7. To the centrifuge tube, add 50 mL of methanol. Swirl well or vortex briefly to loosen the compacted soil. Repeat steps 4-6, combining the extracts in the volumetric flask.
8. Rinse the filter and glassware appropriately with methanol. Add 65.0 mL of NanoPure water to the volumetric flask and adjust the flask to the 200-mL line using methanol.
9. Prepare a Bakerbond 1 gram C<sub>18</sub> SPE column by washing with 2-3 column volumes of methanol followed by 2-3 column volumes of NanoPure water. Do not allow the column to elute to dryness.
10. Volumetrically transfer the requisite volume of each extract to its respective column.
11. Drain the solvent to just above the surface of the C<sub>18</sub> packing material. Wash the column with three 5 mL portions of 70% CH<sub>3</sub>OH : 30% H<sub>2</sub>O. Discard the eluate. Allow the column to dry under vacuum for approximately 15 minutes.
12. To the column, add 5 mL of methanol and elute the residues into a 15-mL graduated centrifuge tube.
13. Place the centrifuge tube in a water bath maintained at 40°C. Evaporate the solvent to dryness under a gentle stream of nitrogen.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 35 -

- 14 Reconstitute the milbemectin residues in the requisite volume of 1% isopropanol (IPA)-99% hexane Vortex each tube for approximately 30 seconds.
- 15
  - A. For analysis of 8,9Z M.A<sub>3</sub>, 8,9Z M.A<sub>4</sub>, transfer the reconstituted extract to an autosampler vial and fractionate using normal phase HPLC. Following evaporation (under nitrogen) of the collected fractions, derivatize the residues according to steps 16-18. See Figure 2.
  - B. For analysis of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub>, volumetrically transfer an aliquot of the reconstituted solution (A) to a second centrifuge tube (B). Evaporate the aliquot (B) to dryness and proceed with step 16.
  - C. For analysis of 27-keto M.A<sub>3</sub> and 27-keto M.A<sub>4</sub>, evaporate (under nitrogen) the remaining solution (A) to dryness and reconstitute to volume with 1% isopropanol (IPA)-99% hexane Vortex each tube for approximately 30 seconds. Analyze by normal phase HPLC.
16. To the dry centrifuge tubes (B), add 1 mL of 0.5M triethylamine (TEA) in benzene and 100 µL of trifluoroacetic anhydride Vortex each tube for a few seconds. Place the tubes in a waterbath (maintained at 40°C) for thirty minutes.
17. Remove the tubes from the waterbath and add 50 µL of 100% TEA. Vortex the solution for a few seconds. Return the tubes to the waterbath and evaporate the solutions under a gentle stream of nitrogen. The solutions will not evaporate to dryness; the final volume achieved should be approximately 200 µL.
18. To each sample, adjust the final volume (using the graduations on the centrifuge tube an/or transferring to volumetric flasks) with 95% CH<sub>3</sub>CN: 5% H<sub>2</sub>O. Vortex the solutions and transfer the final extracts to autosampler vials for analysis.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 36-

Figure 2 Analytical method flow chart for the fractionation of the parent and photoproducts of milbemectin

**METHOD OUTLINE FOR THE FRACTIONATION OF  
M.A<sub>3</sub>, M.A<sub>4</sub>, 8,9Z-M.A<sub>3</sub> AND 8,9Z-M.A<sub>4</sub> IN FIELD SOILS**

1. Follow method outline through #14. (Figure 1)
2. Transfer reconstituted extract to an auto sampler vial and fractionate using normal phase HPLC.
3. Fractionation: set up HPLC. Test inject a standard 2 times to find the correct retention times.
  - a. Once the retention times are established, start the sequence.
  - b. For each sample, collect manually in individual labeled graduated centrifuge tubes (15 mL) the photoproducts (earlier retention time) and the parent material (later retention time). i.e.; If the sample set consists of 10 samples then 20 centrifuge tubes will be needed.
4. Evaporate all fractions under nitrogen and derivatize according to the method outline (Steps 16 through 18, Figure 1.)

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 37 -

**Flow Chart for Processing of Milbemectin and Potential  
Metabolites in Soil**

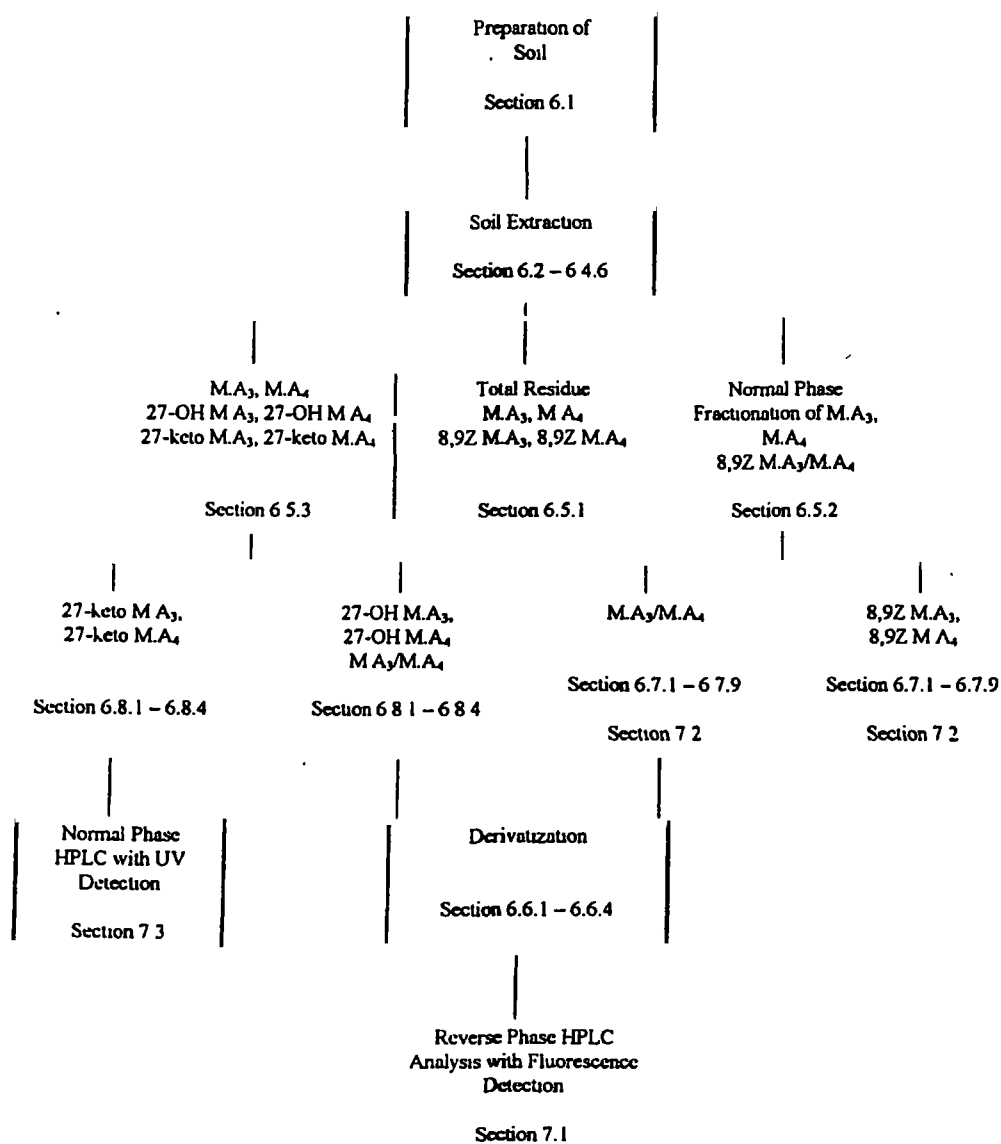


Figure 3 Flow chart for Processing of Milbemectin and Potential Metabolites in Soil

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 38 -

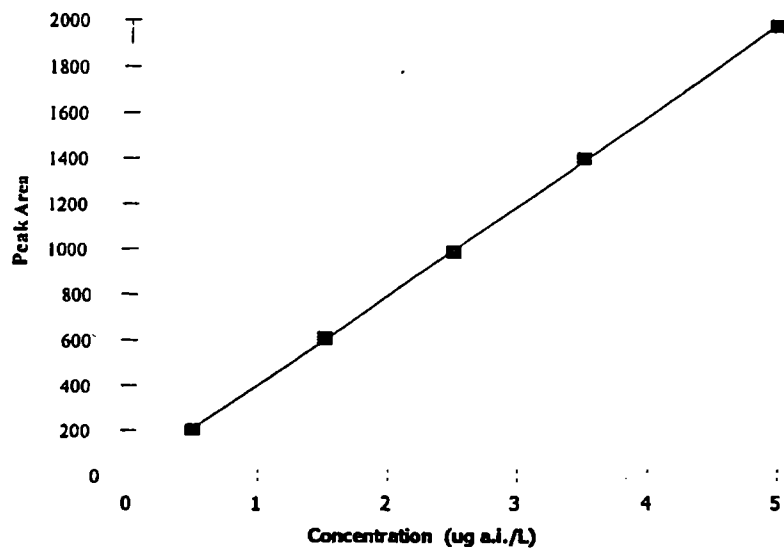


Figure 4 A representative calibration curve for M.A<sub>1</sub>. Slope=392.7285; Intercept=11.93564;  
 $r^2=0.9999$

WILDLIFE INTERNATIONAL LTD

PROJECT NO 420C-106

- 39 -

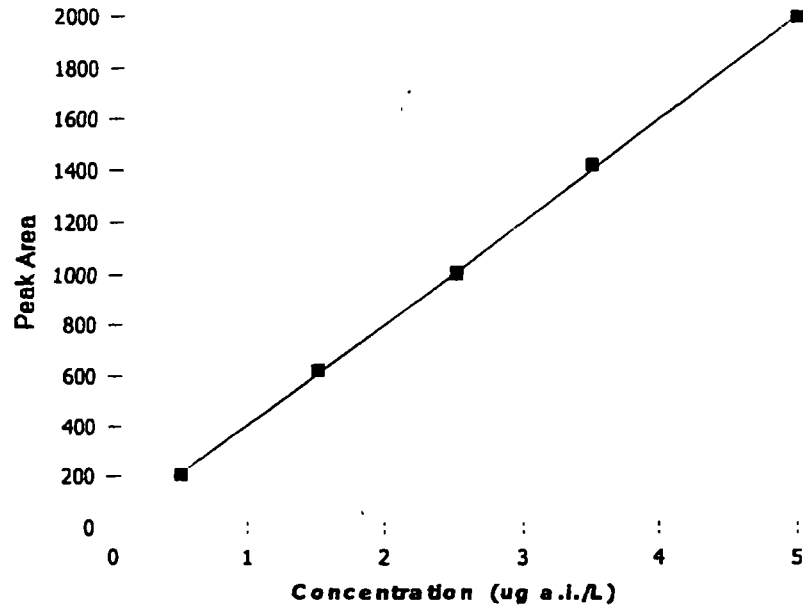


Figure 5. A representative calibration curve for M.A. Slope=397.9923; Intercept=12.92327;  
 $r^2=0.9997$



WILDLIFE INTERNATIONAL LTD

PROJECT NO.: 420C-106

-40-

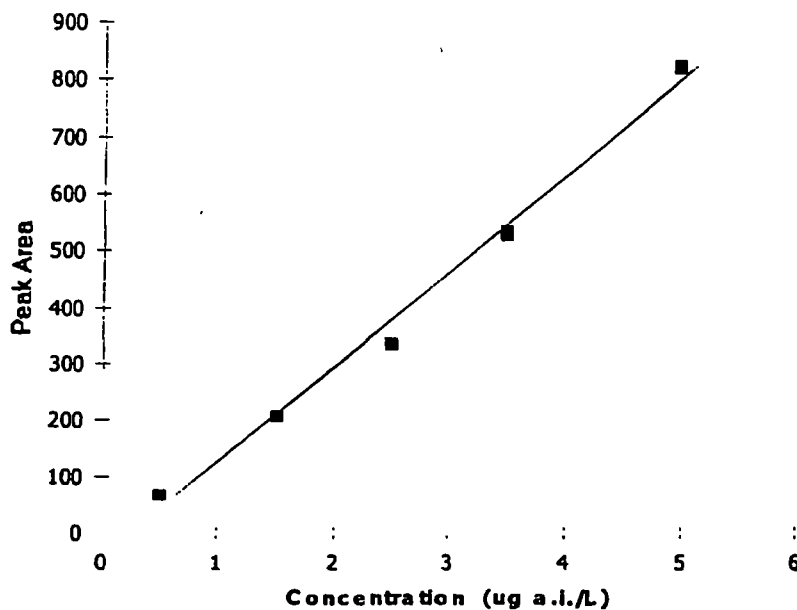


Figure 6. A representative calibration curve for 27-Hydroxy M.A., Slope=167.6254; Intercept=-40 32770;  $r^2=0.9925$

WILDLIFE INTERNATIONAL LTD

PROJECT NO.: 420C-106

-41-

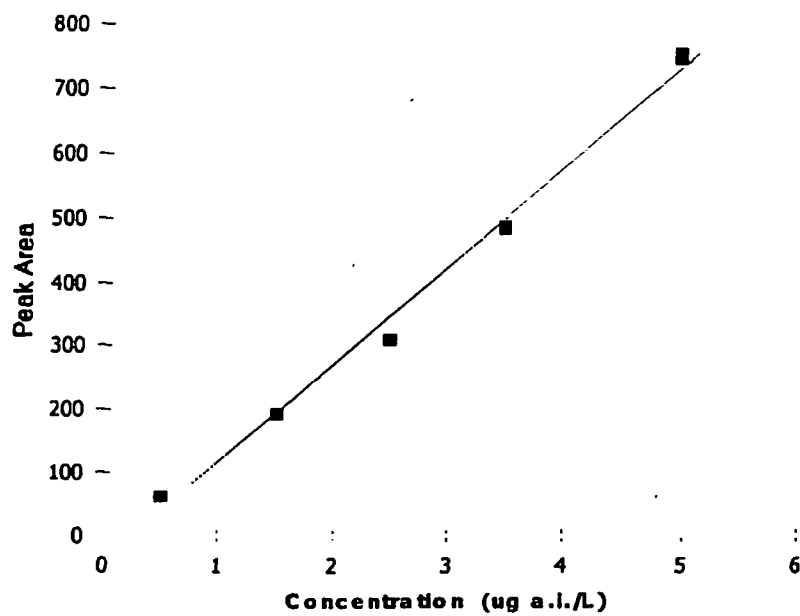


Figure 7. A representative calibration curve for 27-Hydroxy M.A. Slope=152.6734; Intercept=-33.91652;  $r^2=0.9928$

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 420C-106

-42-

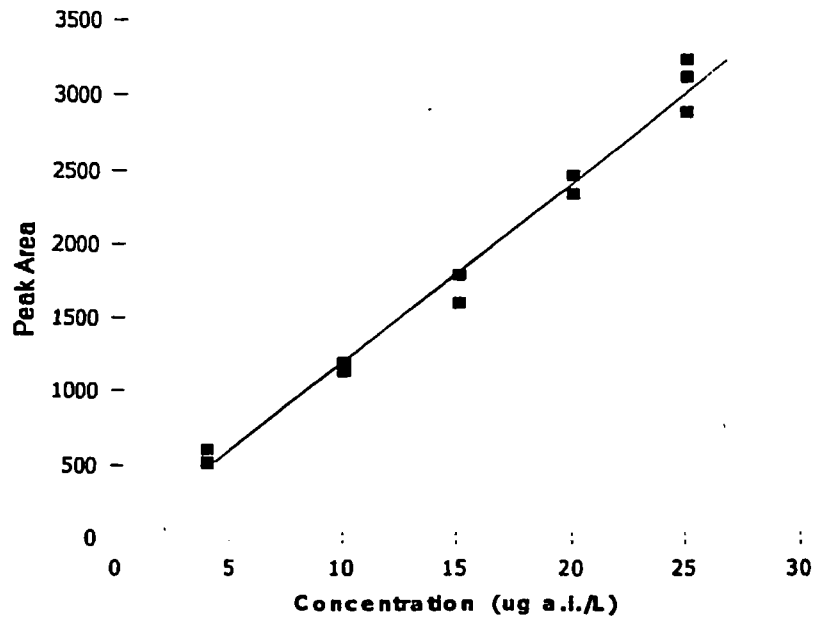


Figure 8. A representative calibration curve for 27-Keto M A<sub>3</sub>. Slope=120.9012; Intercept=-2.98079;  $r^2=0.9868$

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 43 -

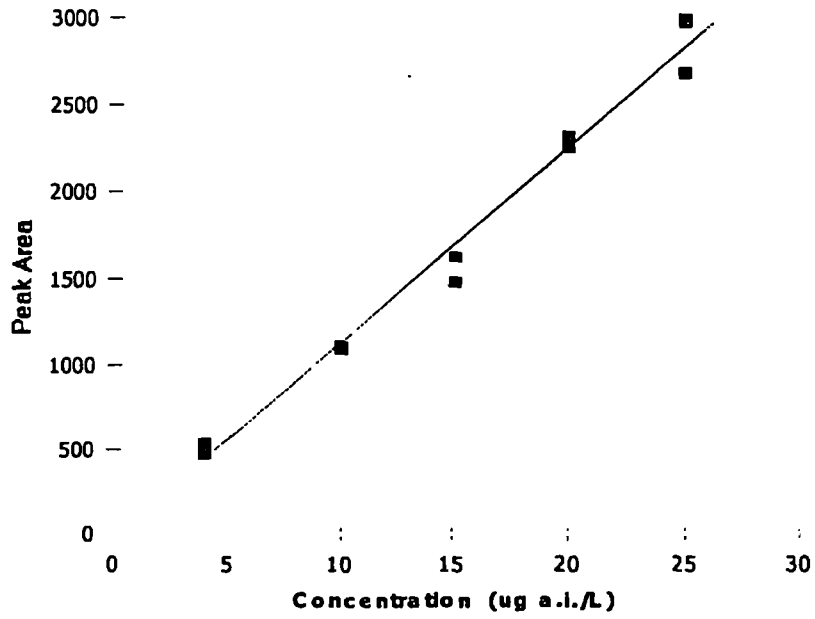


Figure 9 A representative calibration curve for 27-Keto M.A. Slope=113.9810, Intercept=-11.61044;  $r^2=0.9860$

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 44 -

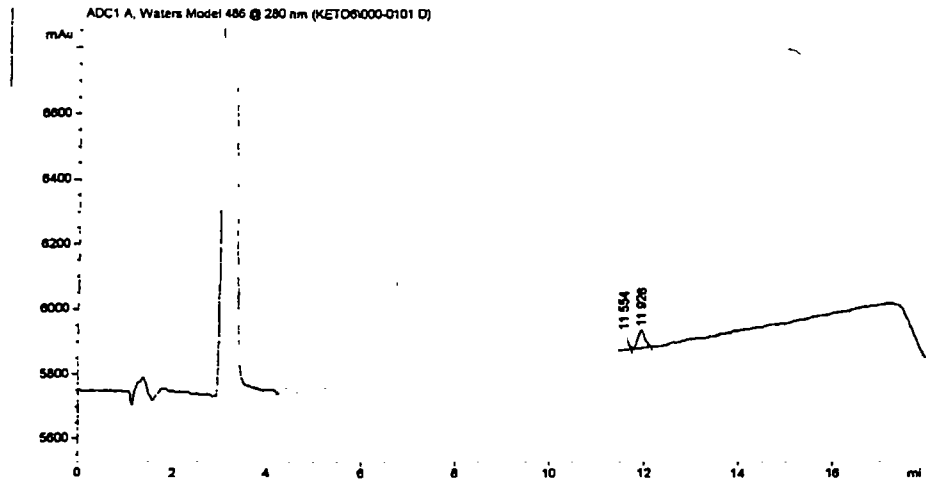


Figure 10 A representative chromatogram of a 5.00 µg a.i./L 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub> combined standard.

WILDLIFE INTERNATIONAL LTD

PROJECT NO 420C-106

- 45 -

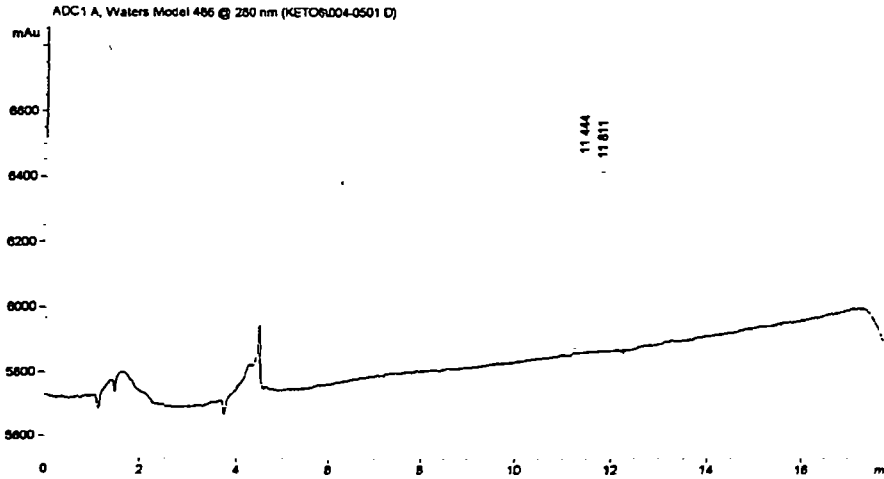


Figure 11. A representative chromatogram of a 50.0 µg a.i./L 27-Keto M.A<sub>3</sub> and 27-Keto M.A<sub>4</sub> combined standard

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 46 -

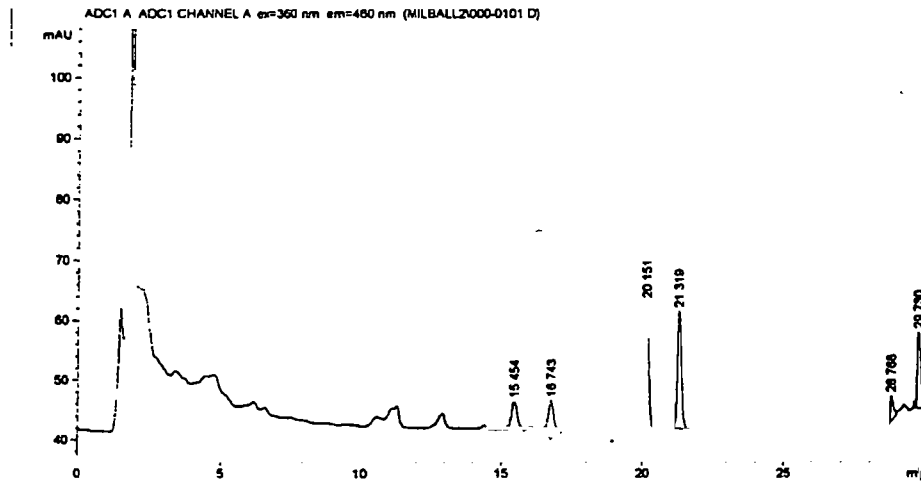


Figure 12 A representative chromatogram of a 0.50  $\mu\text{g}$  a.i./L M.A<sub>3</sub>, 27-Hydroxy M.A<sub>3</sub>, M.A<sub>4</sub> and 27-Hydroxy M.A<sub>4</sub> combined standard

WILDLIFE INTERNATIONAL LTD

PROJECT NO 420C-106

- 47 -

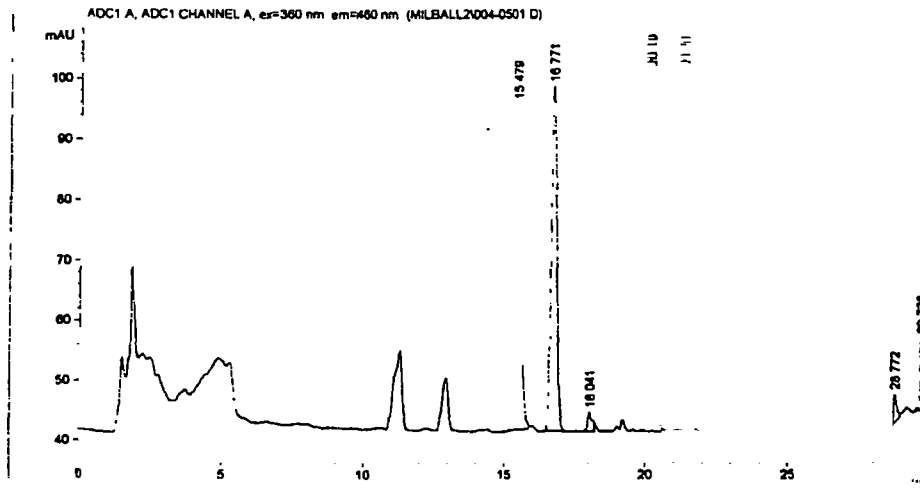


Figure 13. A representative chromatogram of a 5.00 µg a.i./L M.A<sub>3</sub>, 27-Hydroxy M.A<sub>3</sub>, M.A<sub>4</sub> and 27-Hydroxy M.A<sub>4</sub> combined standard



WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 48 -

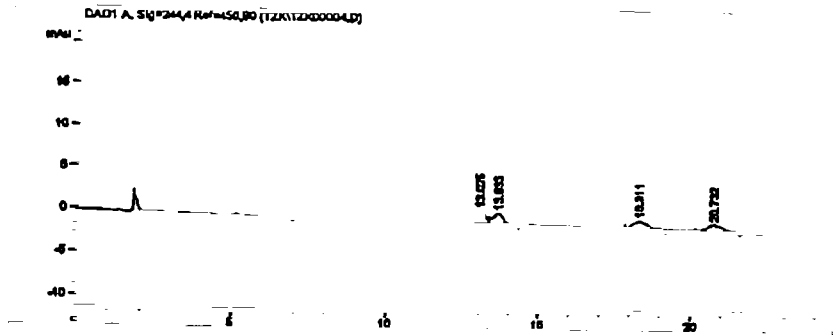


Figure 14. A representative normal phase chromatogram of the fractionation of 8,9Z-M.A<sub>3</sub> and 8,9Z-M.A<sub>4</sub> from M.A<sub>3</sub> and M.A<sub>4</sub>

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 49 -

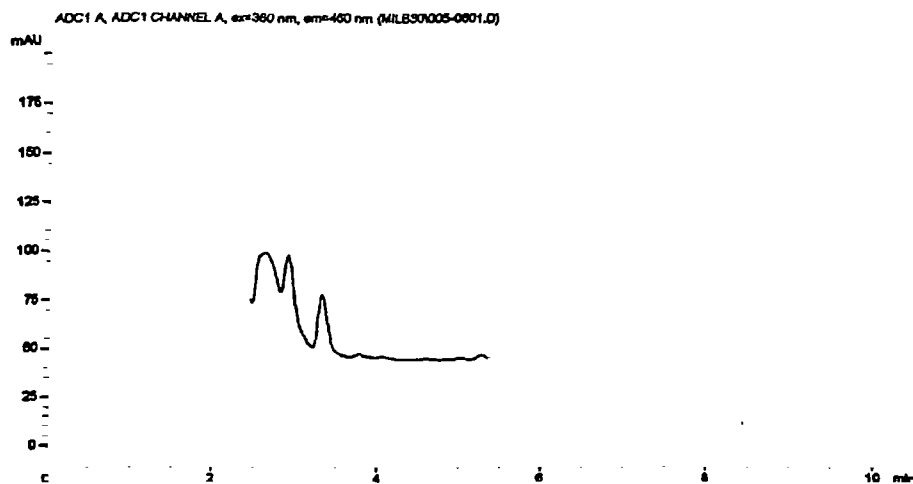


Figure 15. A representative chromatogram of the  $M.A_3/M.A_4$  analysis of the reagent blank from the fractionation of 8,9Z- $M.A_3$ /8,9Z- $M.A_4$  and  $M.A_3/M.A_4$  from New York sandy loam.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 50 -

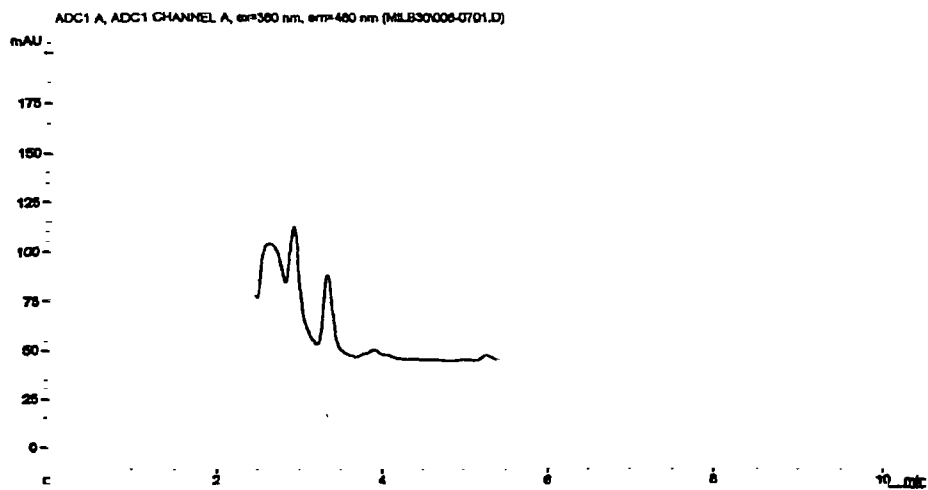


Figure 16. A representative chromatogram of the  $M.A_3/M.A_4$  analysis of the matrix blank from the fractionation of 8,9Z- $M.A_3$ /8,9Z- $M.A_4$  and  $M.A_3/M.A_4$  from New York sandy loam.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 51 -

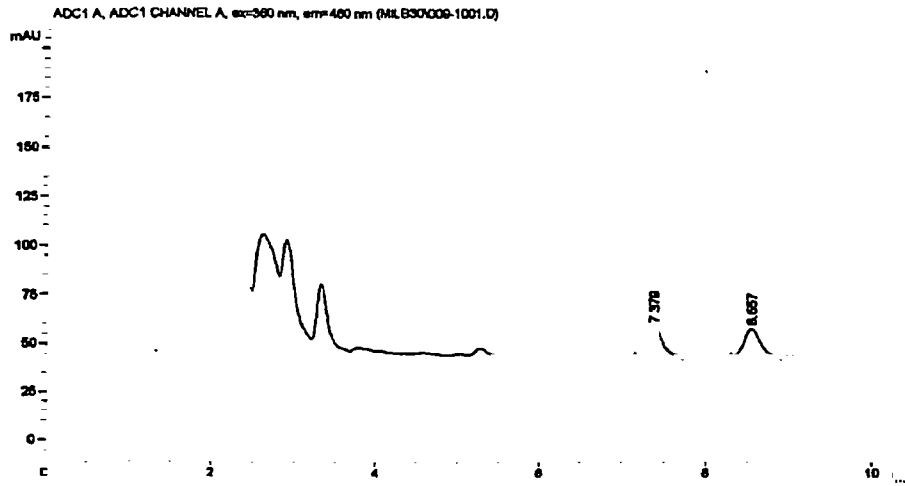


Figure 17. A representative chromatogram of the M.A<sub>3</sub>/M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 2.00 µg a.i./L with 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub> and M.A<sub>3</sub>/M.A<sub>4</sub>.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 52 -

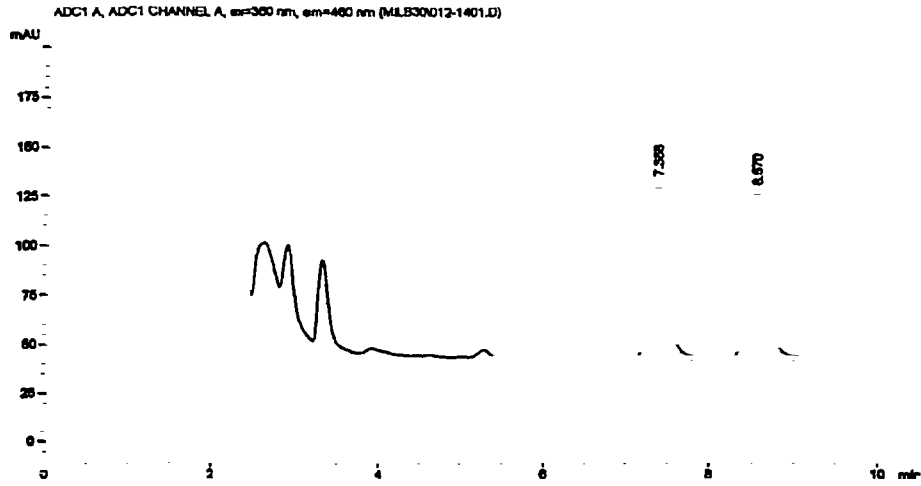


Figure 18. A representative chromatogram of the  $M.A_3/M.A_4$  fractionation analysis of New York sandy loam fortified at  $10.0 \mu\text{g a.i./L}$  with 8,9Z- $M.A_3$ /8,9Z- $M.A_4$  and  $M.A_3/M.A_4$ .

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 53 -

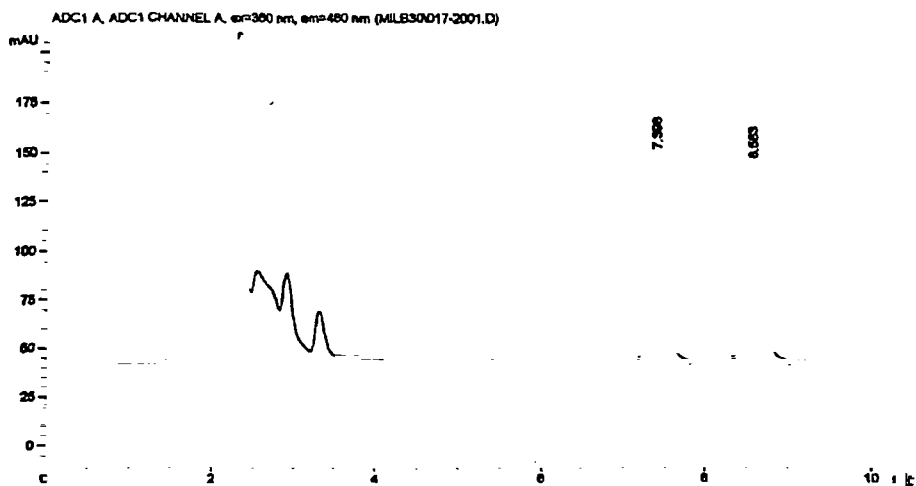


Figure 19. A representative chromatogram of the M.A<sub>3</sub>/M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 50.0  $\mu\text{g a.i./L}$  with 8,9Z-M.A<sub>3</sub>/8,9Z-M.A<sub>4</sub> and M.A<sub>3</sub>/M.A<sub>4</sub>.

WILDLIFE INTERNATIONAL LTD

PROJECT NO. 420C-106

- 54 -

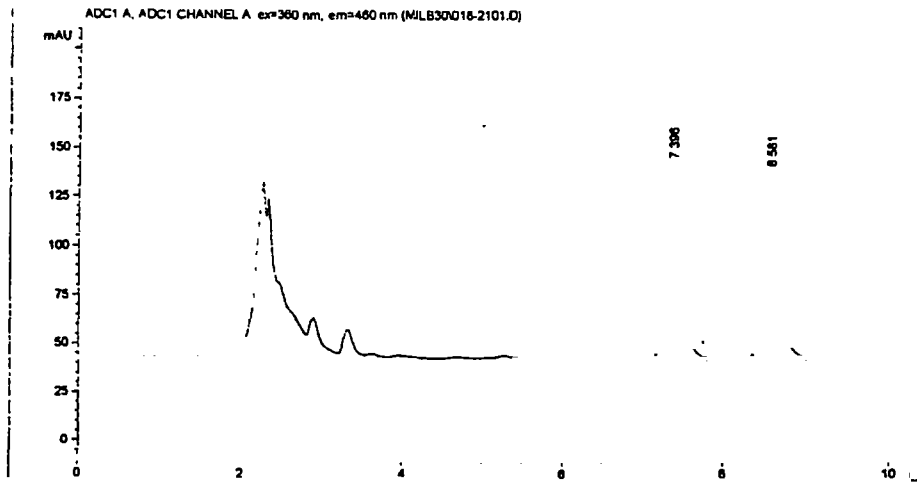


Figure 20. A representative chromatogram of the M.A<sub>7</sub>/M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with 8,9Z-M.A<sub>7</sub>/8,9Z-M.A<sub>4</sub> and M.A<sub>7</sub>/M.A<sub>4</sub>.

WILDLIFE INTERNATIONAL LTD

PROJECT NO: 420C-106

- 55 -

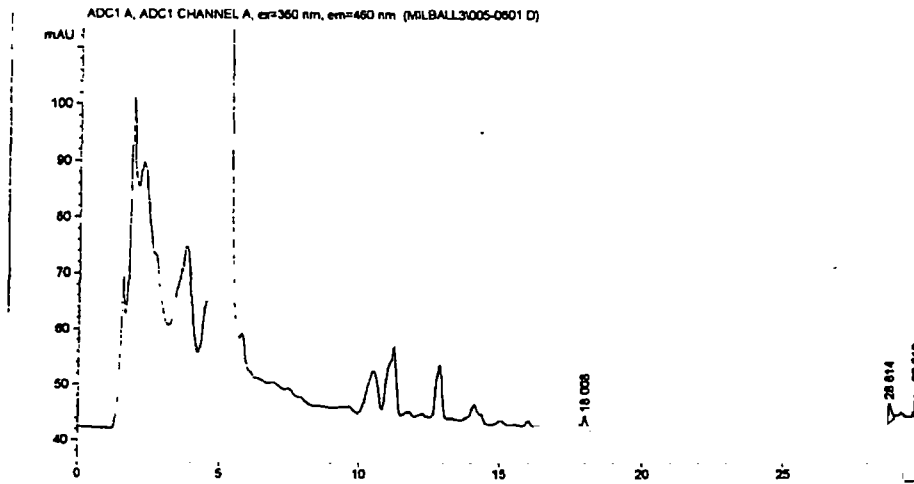


Figure 21. A representative chromatogram of the M.A<sub>3</sub>/M.A<sub>4</sub> and 27-Hydroxy M.A<sub>3</sub>/27-Hydroxy M.A<sub>4</sub> analysis of the reagent blank from the fractionation of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> from New York sandy loam.



WILDLIFE INTERNATIONAL LTD

PROJECT NO. 420C-106

- 56 -

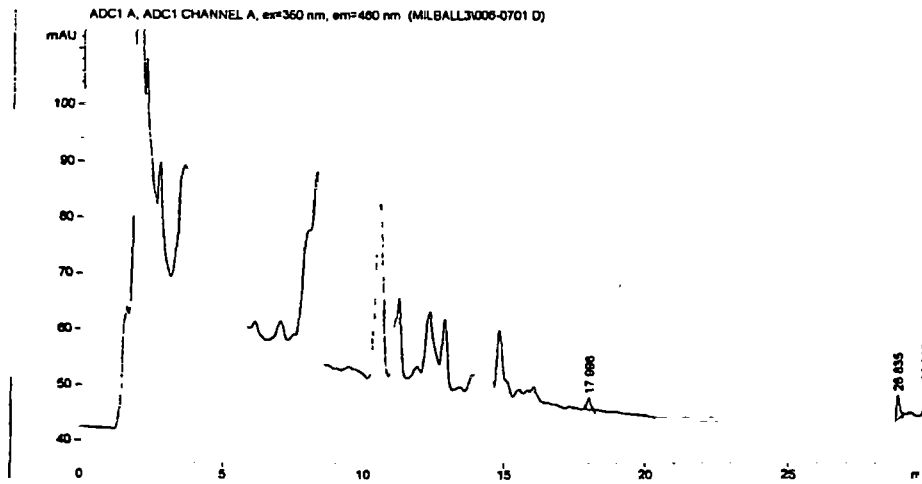


Figure 22. A representative chromatogram of the M.A<sub>3</sub>/M.A<sub>4</sub> and 27-Hydroxy M.A<sub>3</sub>/27-Hydroxy M.A<sub>4</sub> analysis of the matrix blank from the fractionation of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> from New York sandy loam.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO · 420C-106

- 57 -

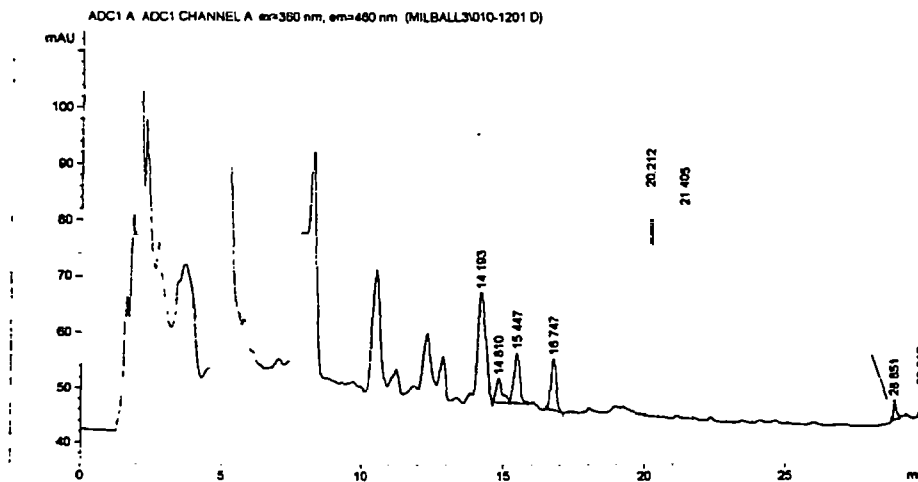


Figure 23. A representative chromatogram of the M.A<sub>3</sub>/M.A<sub>4</sub> and 27-Hydroxy M.A<sub>3</sub>/27-Hydroxy M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 2.00 µg a.i./L with M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub>

WILDLIFE INTERNATIONAL LTD.

PROJECT NO 420C-106

- 58 -

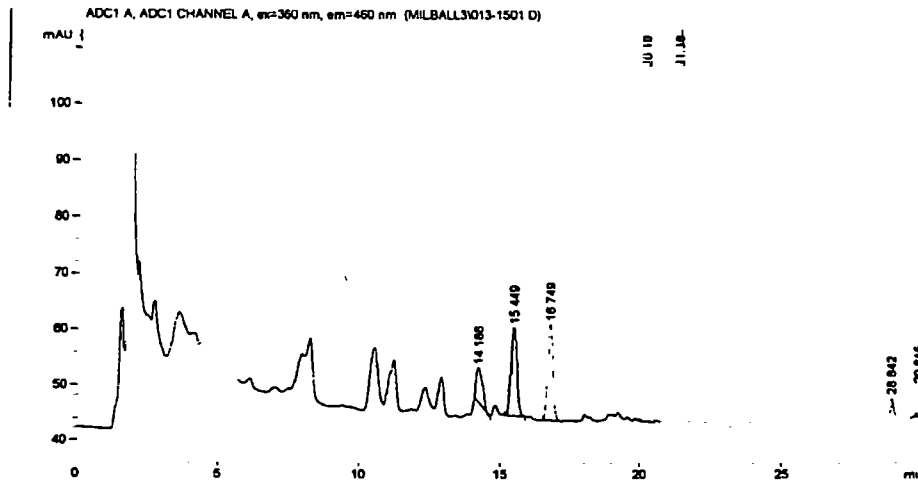


Figure 24. A representative chromatogram of the M<sub>A3</sub>/M<sub>A4</sub> and 27-Hydroxy M<sub>A3</sub>/27-Hydroxy M<sub>A4</sub> fractionation analysis of New York sandy loam fortified at 10.0 µg a.i./L with M<sub>A3</sub>/M<sub>A4</sub>, 27-hydroxy M<sub>A3</sub>/27-hydroxy M<sub>A4</sub> and 27-keto M<sub>A3</sub>/27-keto M<sub>A4</sub>

WILDLIFE INTERNATIONAL LTD.

PROJECT NO : 020C-106

- 59 -

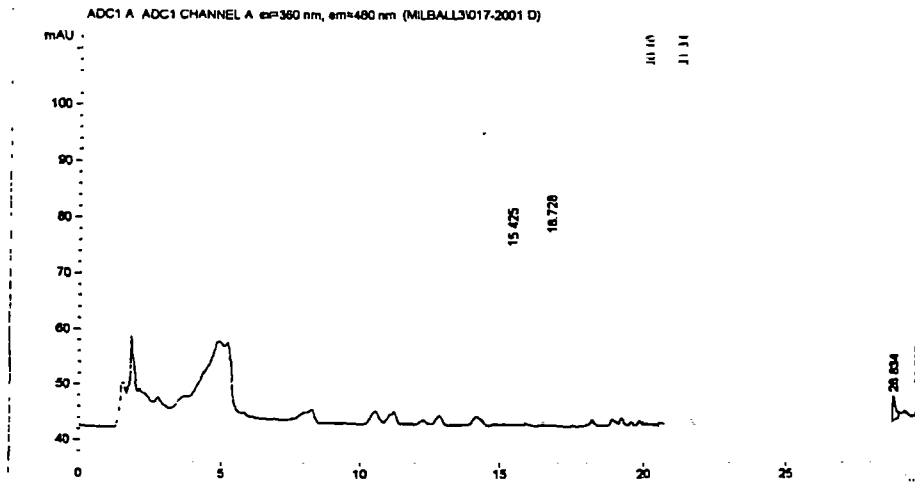


Figure 25. A representative chromatogram of the M<sub>3</sub>/M<sub>4</sub> and 27-Hydroxy M<sub>3</sub>/27-Hydroxy M<sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with M<sub>3</sub>/M<sub>4</sub>, 27-hydroxy M<sub>3</sub>/27-hydroxy M<sub>4</sub> and 27-keto M<sub>3</sub>/27-keto M<sub>4</sub>

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 60 -

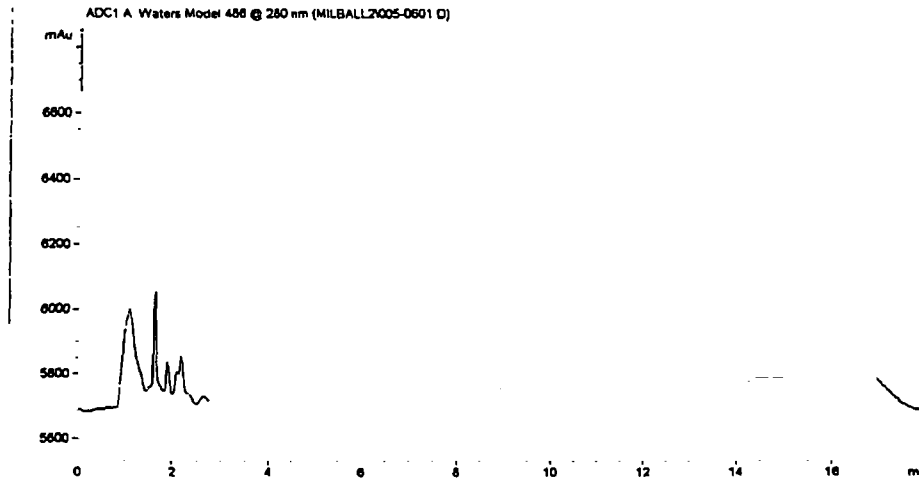


Figure 26 A representative chromatogram of the 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> analysis of the reagent blank from the fractionation of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> from New York sandy loam.

WILDLIFE INTERNATIONAL LTD

PROJECT NO · 420C-106

- 61 -

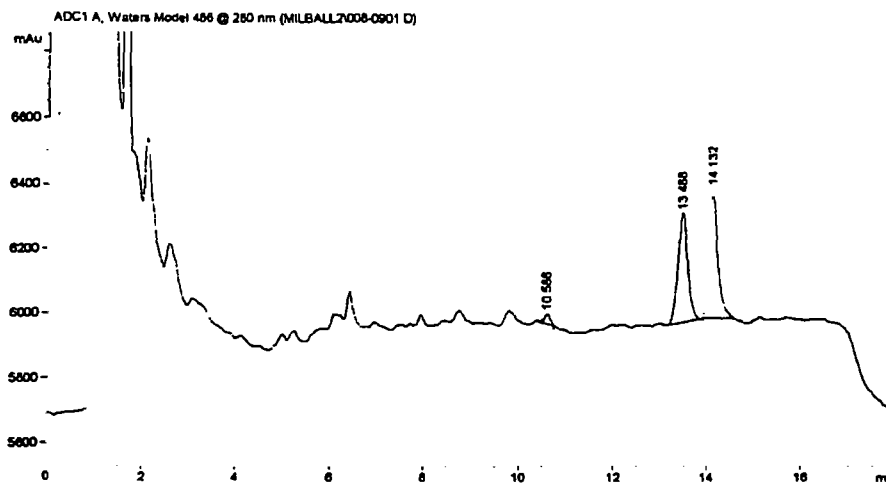


Figure 27. A representative chromatogram of the 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> analysis of the matrix blank from the fractionation of M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> from New York sandy loam

WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 62 -

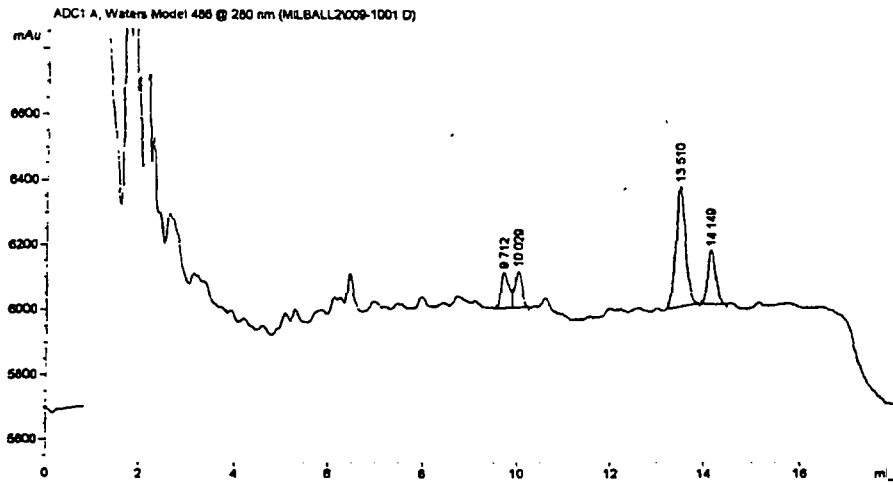


Figure 28. A representative chromatogram of the 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 5.00 µg a.i./L with M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub>.

WILDLIFE INTERNATIONAL LTD

PROJECT NO · 420C-106

- 63 -

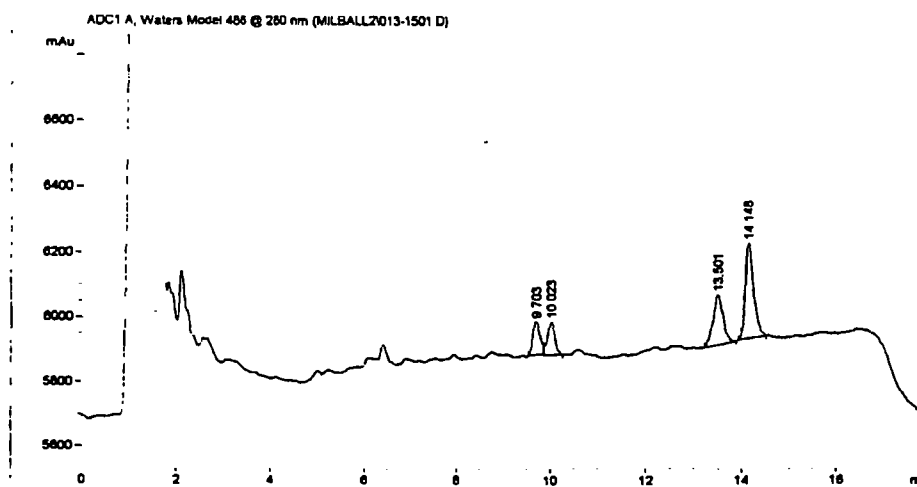


Figure 29. A representative chromatogram of the 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 10.0µg a i./L with M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub>.



WILDLIFE INTERNATIONAL LTD.

PROJECT NO.: 420C-106

- 64 -

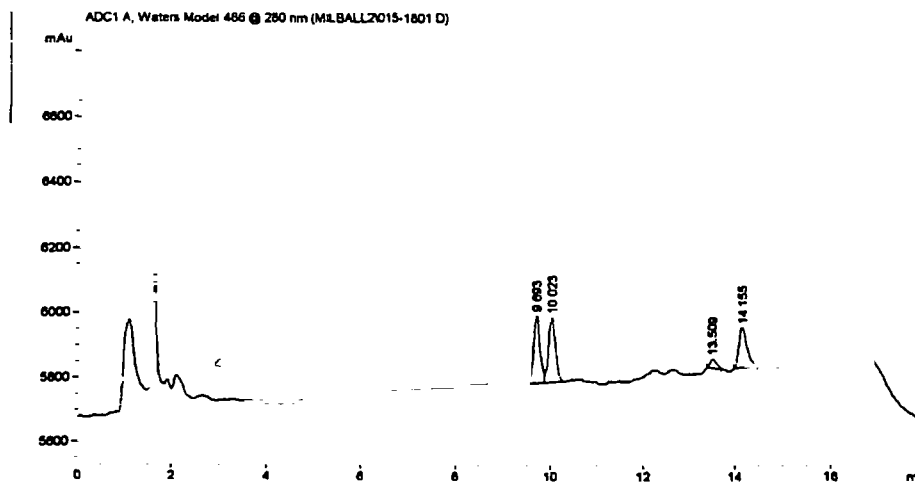


Figure 30. A representative chromatogram of the 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub> fractionation analysis of New York sandy loam fortified at 100 µg a.i./L with M.A<sub>3</sub>/M.A<sub>4</sub>, 27-hydroxy M.A<sub>3</sub>/27-hydroxy M.A<sub>4</sub> and 27-keto M.A<sub>3</sub>/27-keto M.A<sub>4</sub>.

WILDLIFE INTERNATIONAL LTD.

PROJECT NO . 420C-106

- 65 -

## Appendix 1

## Soil Characterization for Milbemectin

	NY
%Sand	65.6
%Silt	24
%Clay	10.4
% Organic Material	3.45
pH	7.0
Cation Exchange Capacity (meq/100g)	7.2
Bulk Density (g/cc)	1.33
Characterization	Sandy Loam

Data acquired from A&amp;L Lab, Ft. Wayne, IN